

Cochlear™

Nucleus® CP910 and CP920 Sound Processor

User Guide



Hear now. And always



This guide is intended for Cochlear implant recipients and their carers using Cochlear™ Nucleus® CP910 or CP920 (“CP900 Series”) Sound Processors. The processor works with your implant to transfer sound to the ear, and is made up of a processing unit, earhook, coil, coil cable and a battery module. The CP920 processor has no accessory socket and is smaller than the CP910 processor.

You can power your processor with a range of battery modules, including a tamper resistant module using disposable batteries, standard and compact rechargeable battery modules, and a LiteWear cable that allows the battery to be worn on the body.

You can control your processor as a stand-alone device using the processor’s buttons, or by using a Cochlear Nucleus CR210 Remote Control or CR230 Remote Assistant. They allow you to control one or two processors at a time and provide extra troubleshooting functions. For more information, please consult your remote user guide.

People with certain types of hearing loss can wear the processor in Hybrid™ mode by adding an acoustic component which sends amplified acoustic sound into the ear canal.

A number of tools and accessories are supplied with your processor.

NOTE

Refer to the Warnings section for warnings and cautions relating to the use of the CP900 Series sound processors, batteries and components.

NOTE

Please also refer to your *Important Information* document for essential advice that applies to Cochlear implant systems.

Symbols used in this guide



Note

Important information or advice.



Tip

Time saving hint.



CAUTION (no harm)

Special care to be taken to ensure safety and effectiveness. Could cause damage to equipment.



WARNING (harmful)

Potential safety hazards and serious adverse reactions. Could cause harm to person.

CP910 Sound Processor

with Standard Rechargeable Battery Module



CP920 Sound Processor

with Compact Rechargeable Battery Module



Contents

Power	5
Batteries	6
Battery life	6
Replace the battery	7
Lock/unlock the battery module to the processor	8
Lock/unlock the tamper resistant battery cover	9
Charge rechargeable battery modules	10
Change disposable batteries	11
Use	13
Turn on and off	14
Change programs	15
Pair with remotes	15
Put on your processor	16
Lock and unlock buttons	18
Wireless accessories	19
Plug-in audio accessories	20
Telecoil	23
Wear	27
Sport and exercise	28
Wear a Mic Lock	29
Wear a Snugfit	30
Wear LiteWear	33
The acoustic component	34
Personalise your processor	40
Travel	41



Store	43
Care	45
Care for your sound processor	46
Care for your battery charger	47
Water, sand and dirt	48
Replace microphone protectors	49
Change earhooks	52
Change the coil cable	57
Change the coil magnet	58
Care for your acoustic component	60
Lights and Beeps	67
Lights	68
Beeps	70
Advanced mode	73
Troubleshoot	76
Warnings	82
For parents and carers	82
Processors and parts	83
Acoustic component	84
Batteries	85
Other information	86

Power

Batteries	6
Battery life	6
Replace the battery	7
Lock/unlock the battery module to the processor	8
Lock/unlock the tamper resistant battery cover	9
Charge rechargeable battery modules	10
Change disposable batteries	11

Batteries

You have a choice of three battery types with the CP900 Series processors:



The Cochlear™ Nucleus® CP900 Series Standard Rechargeable Battery Module.



The Cochlear™ Nucleus® CP900 Series Compact Rechargeable Battery Module.



The Cochlear™ Nucleus® CP900 Series Standard Tamper Resistant Battery Module, which uses two disposable batteries and has a tamper resistant lock.

Battery life


Batteries should be replaced as needed just as you would with any other electronic device. Battery life varies according to the programs used each day, the thickness of skin covering your implant, and the size and type of battery.

The rechargeable battery's lifespan is at least 365 charge cycles. A completely empty rechargeable battery will take approximately four hours to recharge.

Rechargeable batteries may take longer to fully recharge with age. To get the longest life from the rechargeable batteries, always recharge before use.


Replace the battery

Remove the battery

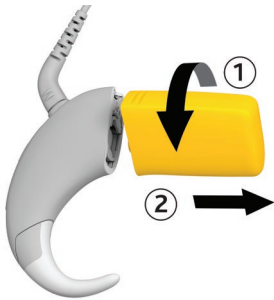

- 1 Twist the battery module as shown to remove it from the processing unit.

Attach the battery

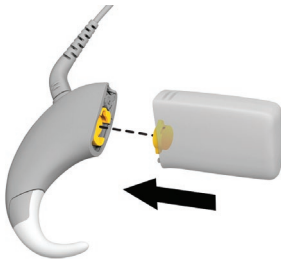

- 1 Align the battery module to the processor socket to fit the parts together.




- 2 Twist the battery module as shown to attach the parts. The processor will turn on automatically.




POWER

Lock/unlock the battery module to the processor


CP910 Processing Unit

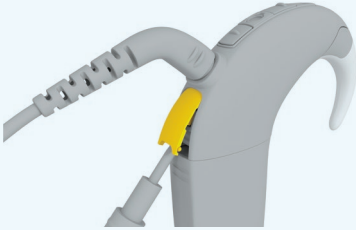


- 1 To increase tamper resistance, carefully lift open the accessory socket cover with your battery cover locking tool. Do not twist the cover.

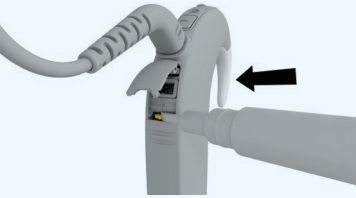

- 2 Push the lock to the far left to lock the battery module.

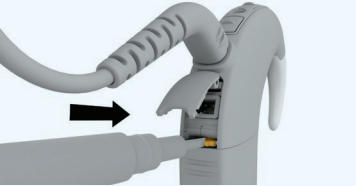

- 3 Push the lock to the far right to unlock the battery module.

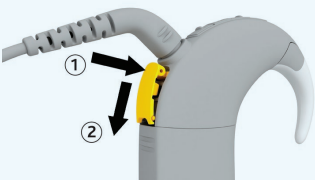



- 4 Carefully close the cover by pushing in the top hinge first.



CAUTION

Always check this lock is unlocked before attaching or removing the battery module.


CP920 Processing Unit

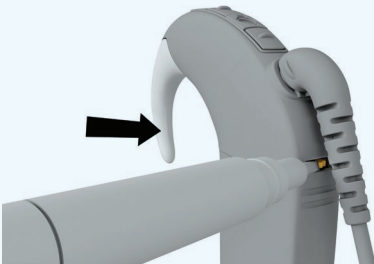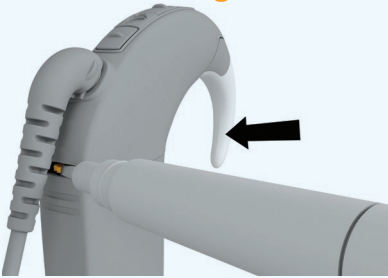

- 1 To increase tamper resistance, push the lock to the far left with your battery cover locking tool to lock the battery module.
- 2 Push the lock to the far right to unlock the battery module.

Lock/unlock the tamper resistant battery cover

The Cochlear™ Nucleus® CP900 Series Tamper Resistant Battery Cover has a tamper resistant lock to help prevent children opening the battery cover.

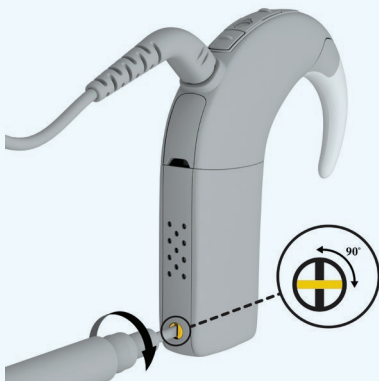

- 1 To lock, turn the locking screw clockwise with the battery cover locking tool until it is in a horizontal position.
- 2 To unlock, turn the locking screw anticlockwise until it is in a vertical position.




Charge rechargeable battery modules







1 Angle and fit the rechargeable battery module into a socket on the Cochlear™ Nucleus® Battery Charger.

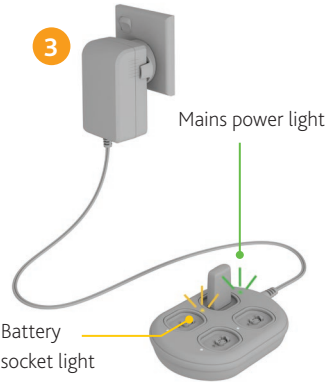

2 Twist clockwise to connect.




3 Connect the battery charger to the power adaptor and plug into mains power.


MAINS POWER LIGHT	WHAT IT MEANS
 Steady green	Battery charger power is on.
Does not light up	Power adaptor is not plugged in or mains power is not available (or, if switched, is not turned on).
BATTERY SOCKET LIGHT	WHAT IT MEANS
 Steady orange	Battery module is recharging.
 Steady green	Battery module is fully charged.
 Flashing orange	There is a problem with the battery module. Try a different socket. Replace the rechargeable battery module.
Does not light up	Battery module is not properly placed, is over-discharged, or there is no power. Check the battery charger is powered on and reconnect the battery charger. If it still does not light up, replace the battery module.

Change disposable batteries


The Cochlear™ Nucleus® CP900 Series Standard Tamper Resistant Battery Module uses two high power zinc air disposable batteries. Cochlear recommends 675 (PR44) zinc air batteries designed for cochlear implant use. Do not use silver oxide or alkaline batteries.

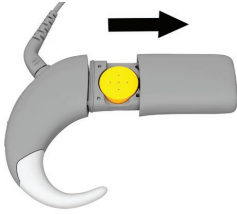


- 1 Pull the unlocked battery cover away from the battery holder.

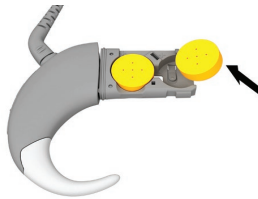


- 2 Remove the batteries from the battery holder.

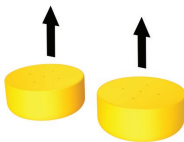

- 3 Remove the new batteries from the packet and let them stand for a few seconds.

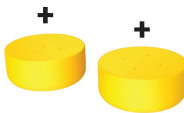

- 4 Insert the batteries into the battery holder with the flat side (positive terminal) facing up.

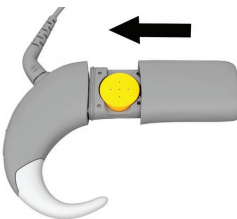

- 5 Replace the battery cover by sliding it up towards the processing unit. The processor will turn on automatically.

See *Batteries* on page 85.



POWER


Use




Turn on and off	14
Change programs	15
Pair with remotes	15
Put on your processor	16
Lock and unlock buttons	18
Wireless accessories	19
Plug-in audio accessories	20
Telecoil	23

Turn on and off

Either connect battery module or short-press lower button to turn on.

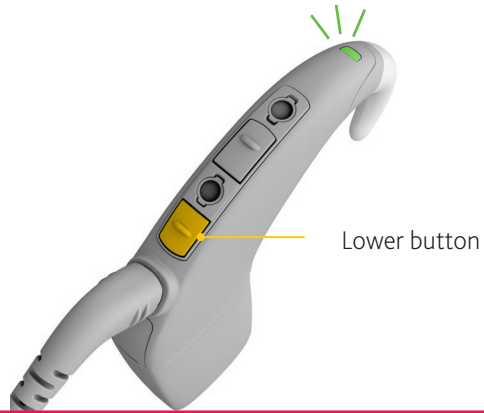Long-press the lower button (about two seconds) to turn off.



INDICATOR LIGHTS	WHAT IT MEANS
 Green flashes	Turning on processor. The number of flashes indicates the number of the current program.
 Steady orange while pressing the lower button	Turning off processor.
 Quick green flashes	Processor flashes while receiving sound from microphones (if set up by your clinician*).

NOTE
Your clinician may set up your processor to turn off automatically when it has been off your implant for more than two minutes.

* See your Cochlear™ Nucleus® CR230 Remote Assistant User Guide for details about modes.

Change programs

You can choose between programs to change the way your sound processor deals with sound e.g. in noisy or quiet places. Usually two programs are all you need, but your clinician can give you up to four programs.



Short-press and release the lower button to switch between programs.

INDICATOR LIGHT

WHAT IT MEANS



Green flashes

Changing the program.
The number of flashes indicates the number of the current program.

Pair with remotes

You need to pair your sound processor to your CR230 Remote Assistant and CR210 Remote Control to use their control and monitoring functions. Please see your remote's User Guide for details.



USE


Put on your processor

1



Place the processor on your ear, letting the coil dangle.

2

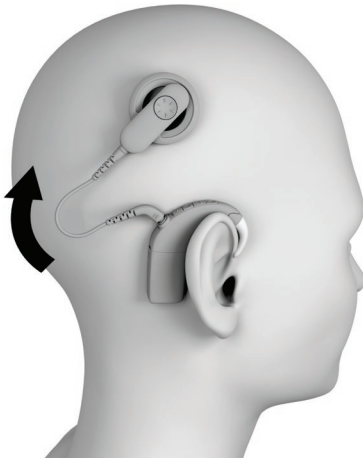

Move the coil sideways and onto your implant.


INDICATOR LIGHTS

WHAT IT MEANS



Flash of orange every second

Processor flashes while coil is off (or connected to the wrong implant).

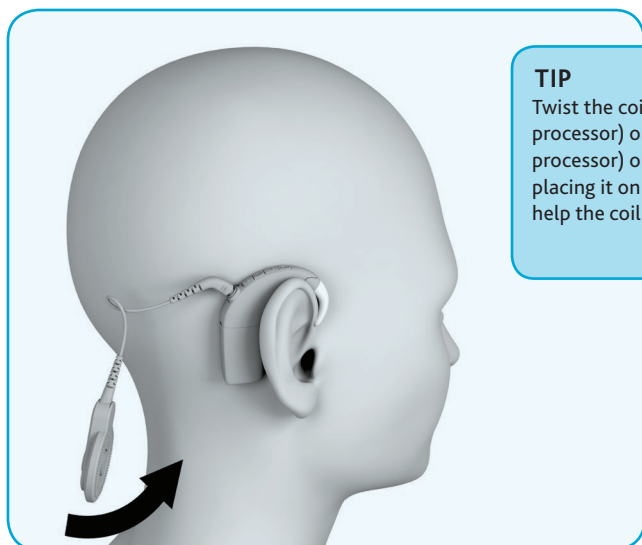

TIP

Twist the coil clockwise (right processor) or anticlockwise (left processor) once **ONLY** before placing it on the head. The twist can help the coil stay in place.






USE

Lock and unlock buttons



Short-press and release both buttons at the same time to lock and unlock the processor.

INDICATOR LIGHTS	WHAT IT MEANS
 Green flash followed by orange	Locking processor buttons.
 Orange flash followed by green	Unlocking processor buttons.
 Orange flash when pressing buttons	Processor buttons are locked.

Wireless accessories

Cochlear True Wireless™ Accessories can wirelessly stream sound to your processor:

- The **Mini Microphone** or **TV Streamer** are controlled from your processor
- You use the **Phone Clip** controls for phone calls.

NOTE

You first need to pair your wireless accessories with your sound processor. See their user guide for details.



USE

Plug-in audio accessories

NOTE

The CP910 sound processor is compatible with the CP800 range of plug-in audio accessories. You can order adaptors to make your CP910 sound processor compatible with Freedom and Euro accessories. Your CP910 sound processor is NOT compatible with the Baha® range of plug-in audio accessories.



The accessory socket on the CP910 processing unit lets you connect a wide range of plug-in audio accessories.



For example, you can choose plug-in audio accessories that help you listen to portable music players, smart phones, computers, TV sets etc. Bilateral recipients can listen in stereo.



You can choose plug-in audio accessories that make conversation easier in noisy environments.




Carers can use a monitor earphone to check that the processor and plug-in audio accessories are receiving sound.

You can order and learn more about accessories at the Cochlear website.

Connect plug-in audio accessories

To connect a plug-in audio accessory to the CP910 processor:

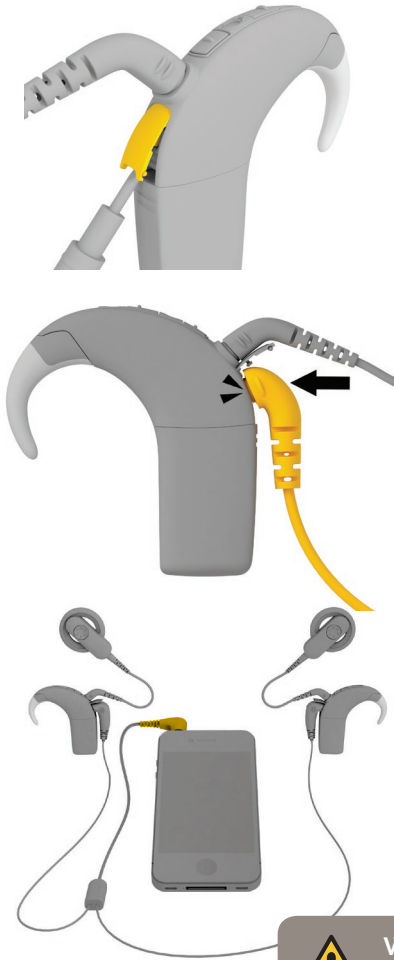

- 1 Carefully lift open the accessory cover with the battery cover locking tool. Do not twist the cover.
- 2 Hold the side grips on the accessory connector and carefully push it into the accessory socket until it clicks into place.
- 3 When you connect a plug-in audio accessory, your processor automatically detects it.
- 4 If applicable, insert the other end of the plug-in audio accessory cable into the compatible device e.g. a portable music player.



WARNING

Always use the mains isolation cable when connecting any personal audio cable (including bilateral) to:

- a mains powered device, e.g. a TV.
- a battery powered device that is connected to mains power (e.g. a laptop connected to an electric outlet for charging).

Do not use excessive force, twist or pull on accessories when connecting them to the processor or at any time while they are connected.






USE

Turn plug-in audio accessories on and off

Plug-in audio accessories turn on automatically when you plug them into your processor.



Short-press the upper button to switch between microphones and the plug-in audio accessory.

INDICATOR LIGHTS	WHAT IT MEANS
 Quick blue flashes	Processor flashes while receiving sound from telecoil/plug-in audio accessory (if set up by your clinician*).
 Long flash of blue	Changing from microphone to telecoil/plug-in audio accessory.
 Long flash of green	Changing from telecoil/plug-in audio accessory to microphones.


* See your Cochlear™ Nucleus® CR230 Remote Assistant User Guide for details about modes.

Telecoil

Listen to hearing loops and hear voices on the telephone with the telecoil. Your clinician can set up your processor with manual telecoil and also auto telecoil. Auto telecoil detects the telecoil signal automatically.

When you use a telecoil-compatible telephone or move in range of a hearing loop, auto telecoil switches to 'receiving'. When there is no signal, e.g. you are out of range of a hearing loop, auto telecoil switches to 'standby' (not receiving).

Whether you choose auto or manual telecoil, you can control the telecoil with your remote. You can turn manual telecoil on and off with your processor.

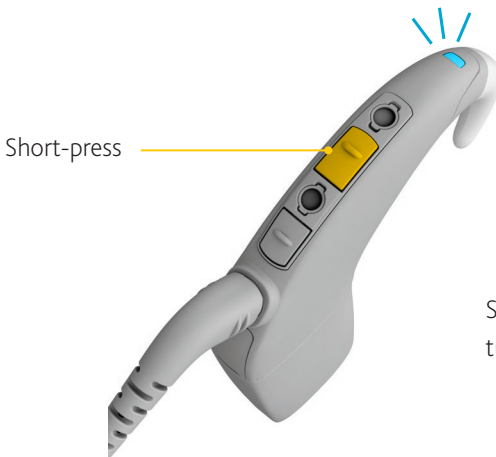

Short-press the upper button to turn manual telecoil on and off.

NOTE

You can also use a Cochlear Wireless Phone Clip to make phone calls. See *Wireless accessories* on page 19.



USE

NOTE

Certain electronic devices or machinery may automatically turn on auto telecoil (you may hear a buzzing noise). If this happens, move away from the device or machinery and wait a few seconds for auto telecoil to turn itself off, or turn it off manually.



TIP

For optimal telecoil performance, position your phone's speaker about 1 cm below the front microphone of your sound processor. It may take a few moments to find the best position.



INDICATOR LIGHTS

WHAT IT MEANS



Quick blue flashes

Processor flashes while receiving sound from telecoil/audio accessory (if set up by your clinician*).



Long flash of blue

Changing from microphone to telecoil/audio accessory.



Long flash of green

Changing from telecoil/audio accessory to microphones.

* See your Cochlear™ Nucleus® CR230 Remote Assistant User Guide for details about modes.



Wear

Sport and exercise	28
Wear a Mic Lock	29
Wear a Snugfit	30
Wear LiteWear	33
The acoustic component	34
Personalise your processor	40
Travel	41


Sport and exercise

TIP

If you want to use your processor while bathing, swimming or showering, ask your clinician about the Cochlear™ Nucleus® Aqua+.



1

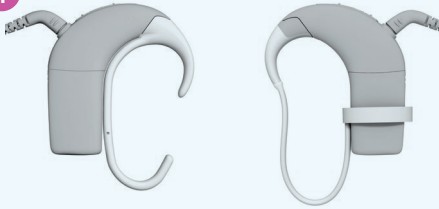

Use Cochlear CP900 Series accessories such as the Snugfit or Mic Lock to help hold your processor in place when you play sport or exercise.

2



After exercise, wipe your processor with a soft cloth to remove sweat or grime.

3



Then check your microphone protectors.



Rechargeable batteries give you the maximum water and sweat resistance.


TIP

Accessories like the Snugfit or Mic Lock can help hold a processor on a child's ear. Tamper resistant accessories such as the tamper resistant earhook, tamper resistant battery module and tamper resistant Snugfit are available. Always ensure these are properly locked together before use.





Wear a Mic Lock

The Cochlear™ Nucleus® Mic Lock helps hold the processor securely on the ear.

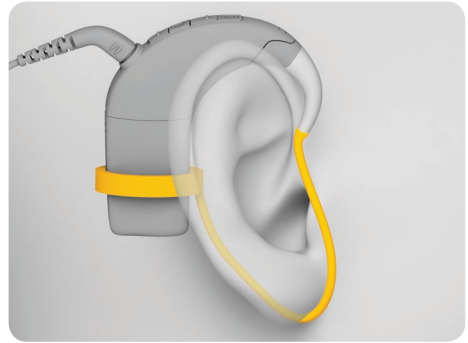

To attach the Mic Lock the first time:

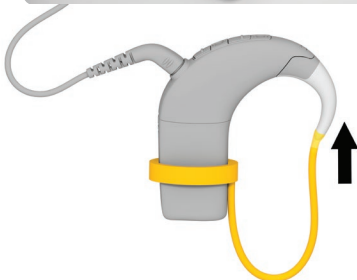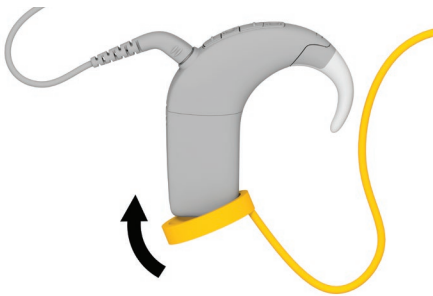

- 1 Push the Mic Lock band onto the processor with the tubing facing forward.
- 2 Place the processor on your ear and bring the tubing up to the earhook.
- 3 Trim the tubing to a length that will hold the processor firmly in place without discomfort. Do not cut it too short.
- 4 Soften the end of the tubing in warm water and push it onto the tip of your earhook.
- 5 Let the tubing cool before wearing.



WEAR

Wear a Snugfit

The Cochlear™ Nucleus® CP900 Series Snugfit holds your sound processor more securely in place than an earhook alone. It gives carers easy access to controls and accessories and is available in large, medium and tamper resistant (small).




TIP

You can wear Snugfits with sound processor covers and coil covers.



Attach a Snugfit

- 1 Remove the standard earhook. Either:
 - Press firmly on the top to click it off,
 - or
 - Twist the earhook off.

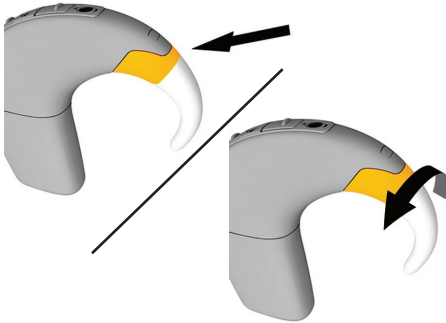

- 2 Click your Snugfit into place.



- 3 Carefully bend the lower part into shape to suit your ear.



Remove a Snugfit


- 1 Remove your Snugfit by clicking it off the processor.




Attach a Tamper Resistant Snugfit




- 1 The smallest Snugfit comes with a Snugfit band to fix it more securely to the processor.

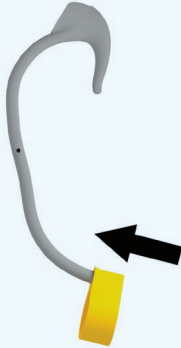

- 2 Thread the Snugfit band onto the Snugfit as shown.

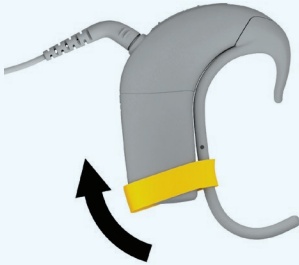

- 3 Click the Snugfit into place on the processor and push the Snugfit band onto the battery module.

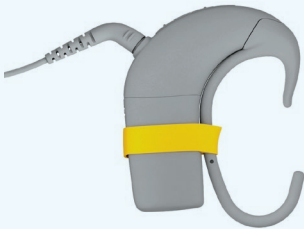

- 4 Push the band securely onto the battery module.

Wear LiteWear

The Cochlear™ Nucleus® CP900 Series LiteWear Cable allows any standard (not compact) battery module to be worn on the body – so only the processing unit is worn on the ear. You can order LiteWear from Cochlear or talk to your clinician.



LiteWear is made up of a LiteWear cable, case and a choice of fixing aids.

NOTE

You can wear a Snugfit with LiteWear to hold the processor more securely on the ear.



TIP


To help prevent tampering, attach the LiteWear case between the shoulder blades of very young children.




The acoustic component

If your sound processor is programmed for use in Hybrid mode, you can wear it with an acoustic component.

Wear your acoustic component

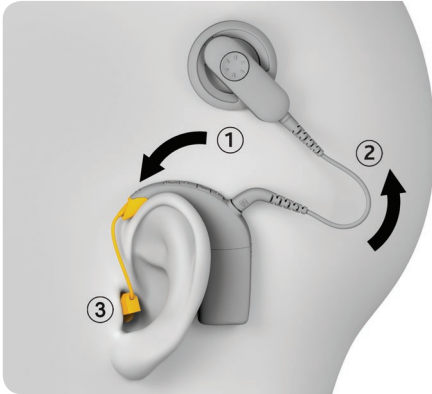

- 1 Hold the coil in one hand and place the processor on your ear with the other.
- 2 Place the coil on your head.
- 3 Place your dome or earmould in your ear as instructed by your clinician.
- 4 When removing your coil and processor, remove your dome or earmould as instructed by your clinician.

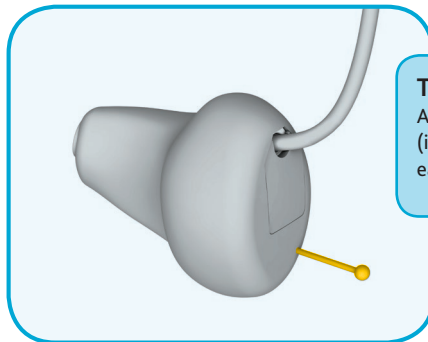


TIP

Always use the removal cord (if present) when taking the earmould out of your ear.



Remove your acoustic component

Wear your processor without the acoustic component by replacing it with a standard earhook. Talk to your clinician before replacing:

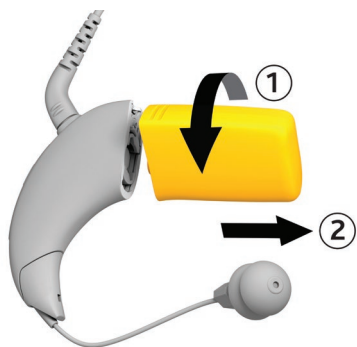

- 1 Remove the battery module from the processing unit by twisting the two parts as shown.



- 2 Locate the opening under the earhook.


See *Acoustic component* on page 84.






CAUTION


Use only the EAC200 Series Screwdriver – not the Cochlear™ Nucleus® Battery Cover Locking Tool – to remove the acoustic component.

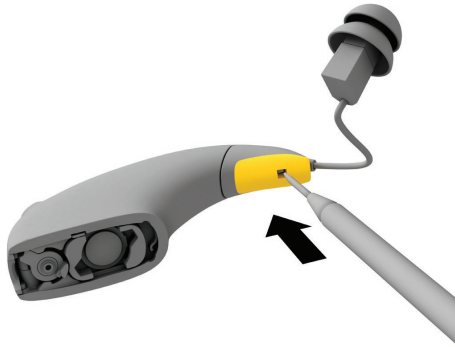


- 3 Insert a Cochlear EAC200 Series Screwdriver into the hole.

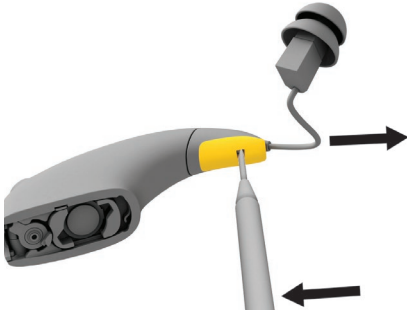

- 4 Carefully push in, and lever the screwdriver towards the processor, until the speaker unit cable disconnects.

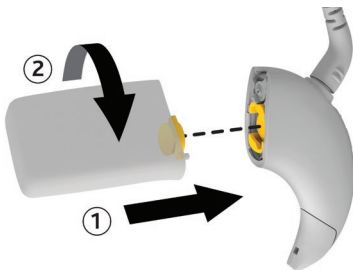


- 5 Re-attach the battery module.

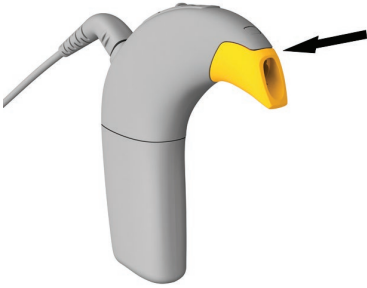

- 6 Press firmly on top of the Cochlear EAC200 Series Earhook with your thumb to click it off the sound processor.



- 7 Insert a new processing unit plug into the speaker unit cable socket. Make sure the plug goes all the way in.

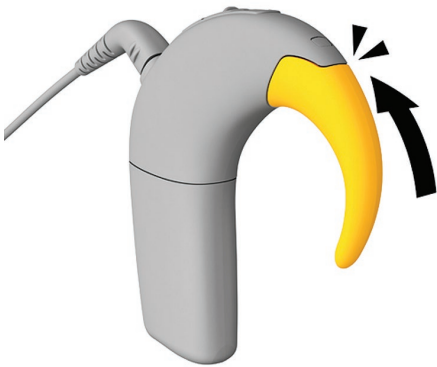


- 8 Click the standard earhook into place.




WEAR

Attach your acoustic component

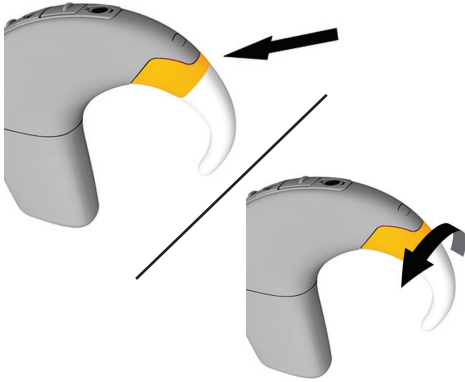

- 1 Remove the standard earhook.
Either:
 - Press firmly on the top to click it off,
 - or
 - Twist the earhook off.

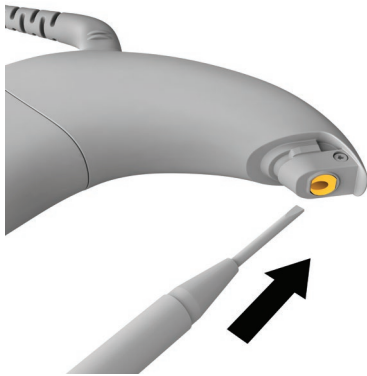


- 2 Insert the Cochlear™ EAC200 Series Screwdriver into the hole in the processor unit plug.




- 3 Use the screwdriver to carefully remove the plug. Do not twist or try to lever the plug out. Discard the plug.



- 4 Click the EAC200 series earhook into place.

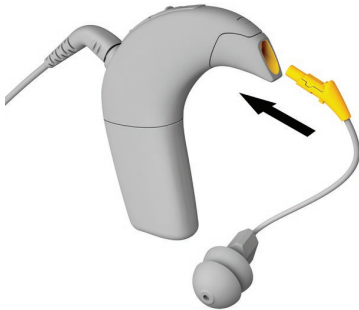

- 5 Align the speaker unit plug with the end of the earhook and carefully slide the plug into the processor.

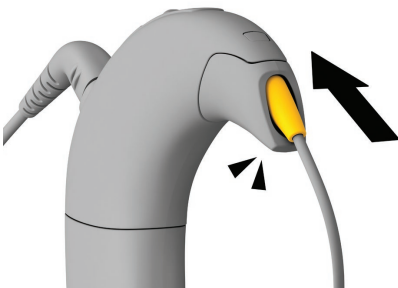

- 6 It should click securely into place. If necessary, carefully push it until you feel the click. Pull gently on the cable to check that it's securely in place.



WEAR

Personalise your processor

You can order Cochlear™ Nucleus® Sound Processor Covers and Coil Covers to personalise your processor and protect it from scratches.

Attach and remove processor covers



- 1 To attach, slide the cover onto the processor.
- 2 To remove, pull the cover away from the processor.

Attach and remove coil covers



- 1 To attach, discard the clear plastic protector on the cover and click the cover into place on the coil.
- 2 To remove, lift the cover off the coil.

Bilateral stickers

If you have two Cochlear implants, you must match the correct processor to the correct implant. Ask your clinician to give you bilateral stickers (red for right, blue for left) to make identifying left and right processors easier.

Travel

- Take a printout from your clinician of your most recent program in case you need help with your processor.
- If you have a back-up sound processor, check that it is programmed correctly and take it with you.
- It's okay to move through metal detectors and full body scanners with your sound processor on. To avoid any possible buzzing sounds in your ear, turn off the telecoil.
- Ask your clinician for a Patient Identification Card. In the unlikely event that your implant sets off a metal detector the ID card will help explain that you have an implanted medical device.
- If you need to remove your sound processor as you move through a security checkpoint, place it in a case in your hand luggage.
- Your sound processor will not interfere with a plane's navigation system so you won't need to turn it off during takeoff and landing. For advice about using your remote control or remote assistant on a plane, refer to your Important Information document.

NOTE

Visit www.cochlear.com/clinic-finder to find the nearest clinic in places you are travelling.




WEAR



Dry aid kit

Store your processor at night in the dry aid kit provided by Cochlear. Store the processor fully assembled for 8 hours for optimal drying effect.

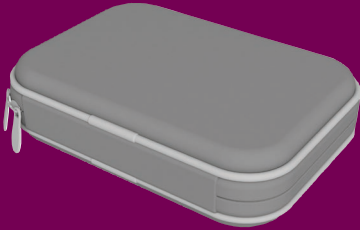

Storage case

For long term storage, remove the batteries and store so they do not touch each other. Storage cases are available from Cochlear.

Store

Care

Care for your sound processor	46
Care for your battery charger	47
Water, sand and dirt	48
Replace microphone protectors	49
Change earhooks	52
Change the coil cable	57
Change the coil magnet	58
Care for your acoustic component	60

Care for your sound processor

Every day

- Check all parts for dirt and moisture. Wipe the processor, coil cable, earhook, covers, Snugfits and Mic Lock cables with a soft dry cloth. (You can leave the earhooks, Snugfits, Mic Locks or coil cables on the processor while cleaning.)
- Keep your processor free from moisture by drying it every night in your dry aid kit.
- Remove the battery module and make sure all the contacts are clean. Carefully tap or blow on them to remove any dirt.
- Make sure the accessory socket and audio accessories are clean. Blow on the accessory socket to remove any dust and clean accessories with a soft dry cloth.
- Check the microphone protectors for signs of dirt or grime and replace if needed. See *Replace microphone protectors* on page 49.

Every month

- Check if earhooks or Snugfits are becoming loose or showing signs of wear. Replace as needed. See *Wear a Snugfit* on page 30, or *Change earhooks* on page 52.

Every two months

- Replace the dry brick in your dry aid kit.

Every three months

- Replace the microphone protectors. See *Replace microphone protectors* on page 49.

NOTE

If your sound processor is subject to increased moisture or dusty conditions, consider cleaning and maintaining parts more often to keep it functioning at its best.



Care for your battery charger

Every day

- Check your battery charger is clean. If you notice any dust or dirt:
 - 1 Unplug the power adaptor and remove any battery modules.
 - 2 Hold the battery charger upside down and tap it gently to remove any dirt from the battery charger sockets. Carefully blowing on the battery charger sockets may also help remove dirt.
 - 3 Wipe the battery charger sockets with a soft dry cloth.
- Cycle batteries in different sockets when recharging to wear the sockets evenly.

If it gets wet

- If the battery charger ever gets splashed with liquid, carefully shake out the liquid and dry the battery charger for 24 hours. Do not use the battery charger until it is dry.



Water, sand and dirt

Your sound processor is resistant to water, sand and dust. However, it is still a precision electronic device so the following precautions should be taken.



If your processor ever gets wet, dry it with a soft cloth, change the microphone protectors and place it in the dry aid kit provided by Cochlear for 8 hours.



If sand or dirt ever enter the processor, shake the components carefully to remove it and change the microphone protectors.



Your processor is protected against failure from dust and temporary immersion in water (IP57 rated) when you wear it with:

- a rechargeable battery module
- a coil and coil cable
- a closed accessory socket
- no acoustic component.



Your processor is protected against failure from dust and splashing water (IP54 rated) when you wear it with:


- a tamper resistant battery module
- a coil and coil cable
- a closed accessory socket
- an acoustic component worn in the ear.

Replace microphone protectors

Replace your microphone protectors every three months, whenever they look dirty or whenever you notice any loss in sound quality. Always replace both microphone protectors at the same time.



Remove microphone protectors

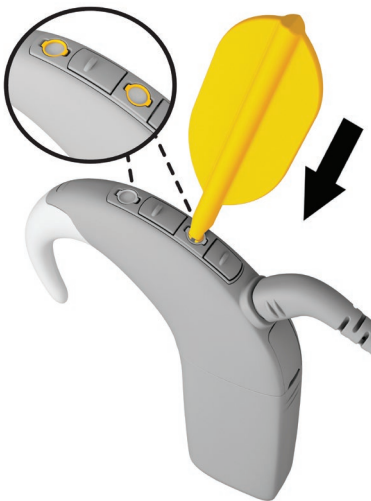

- 1 Use the Cochlear™ Nucleus® Microphone Protector Removal Tool.
- 2 Firmly push the tip of the tool into the middle of the microphone protector.






3 Turn the tool 90° clockwise.




4 Lift out the used microphone protector.

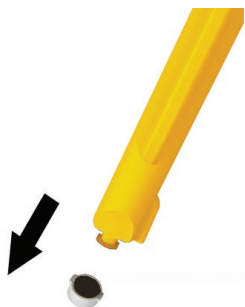

5 Pull the used microphone protector from the tool and discard.

6 Repeat steps 2 to 5 to remove the other microphone protector.

Insert new microphone protectors

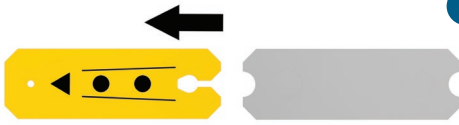

- 1 Pull the Cochlear™ Nucleus® Microphone Protector Applicator out of its sleeve.



- 2 With the microphone protectors facing the processor, gently push the slotted end of the applicator under the coil cable plug.



- 3 Lay the applicator over the processor, aligning the microphone protectors with the microphones.



- 4 Press the microphone protectors all the way into the microphone holes one at a time.




- 5 Remove the applicator carefully, peeling it upwards from the earhook end.




Change earhooks

Cochlear™ Nucleus® CP900 Series Earhooks come in large, medium and small sizes. The smallest earhook can be made tamper resistant.

Change an earhook

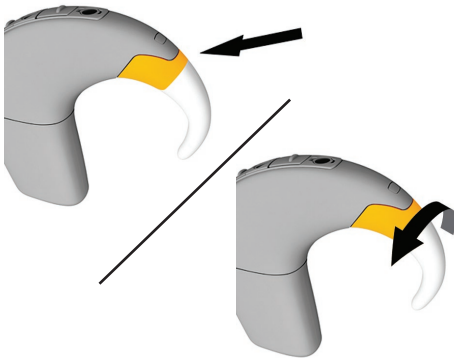

- 1 Remove the standard earhook. Either:
 - Press firmly on the top to click it off,
 - or
 - Twist the earhook off.

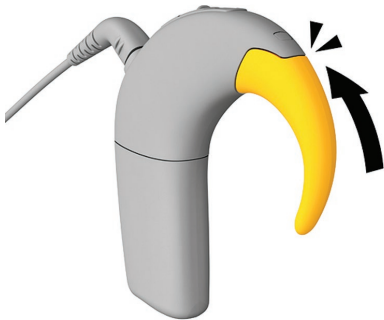

- 2 Click the earhook into place.


NOTE

Remove your earhook only when necessary – it may become loose if it is removed too often.



About the tamper resistant earhook


The tamper resistant earhook can be locked to the processor with a fitting pin to reduce the risk of the earhook becoming a choking hazard for young children. Older children and adults can wear it without the fitting pin.

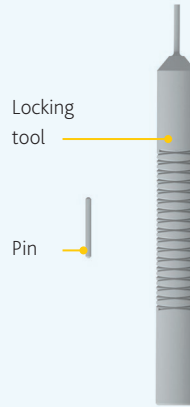


The tamper resistant earhook comes with a locking tool and fitting pins. Use black fitting pins only.




Attach the tamper resistant earhook

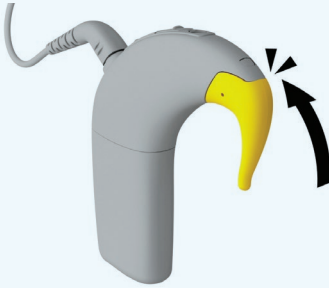


- 1 Click the tamper resistant earhook onto the processor.

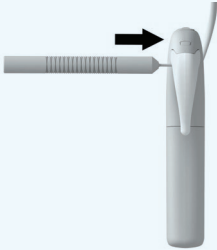


- 2 Insert the pointed end of the locking tool into the hole in the earhook to clear and align the hole.

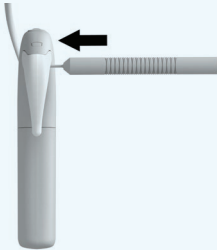

- 3 Flip the processor and insert the pointed end of the locking tool into the hole on the other side of the earhook to clear and align the hole.

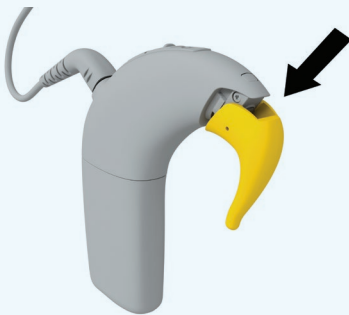

- 4 Press firmly on the top of the earhook with your thumb and click it off the sound processor.

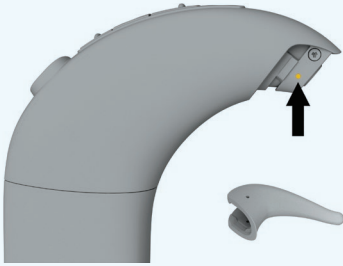


- 5 Make sure the hole is clear. Blow debris away or clear with the pointed end of the tool.


- 6 Replace the earhook.



- 7 Insert a fitting pin into the hole on the blunt end of the locking tool.

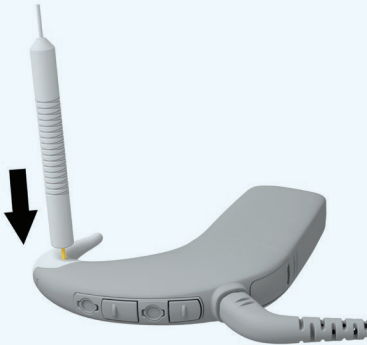


- 8 Hold the processor on a table. Partly insert the fitting pin into the hole on the earhook.

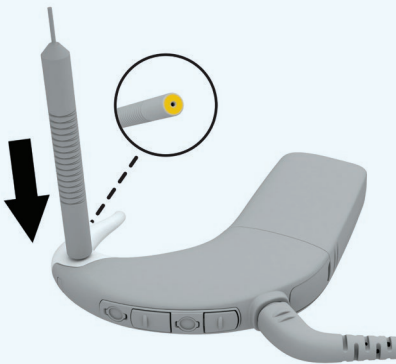

- 9 Use a flat part of the blunt end of the tool to completely push the fitting pin in. The fitting pin should not protrude from either side of the earhook.



Remove the tamper resistant earhook

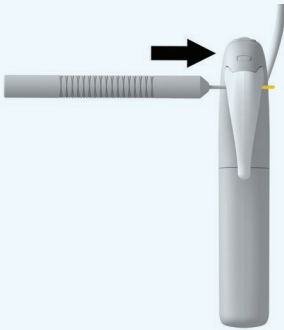


- 1 Insert the pointed end of the locking tool into the hole on the earhook.

- 2 Push the tool into the hole to partly remove the fitting pin. Be gentle as too much force may damage the tool.



- 3 Lift out the earhook fitting pin.

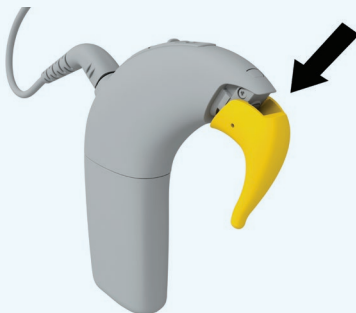

- 4 Press firmly on the top of the earhook with your thumb and click it off the sound processor.

Change the coil cable

The Cochlear™ Nucleus® CP900 Series Coil Cable comes in four lengths so you can choose the most comfortable fit. Only remove the coil cable when you are changing it.




- 1 Hold the coil cable grip and firmly **pull it straight out** of the processor. Do not tug on the flexible part of the coil cable.



CAUTION


Do not twist the coil cable when you pull it out of the processor.

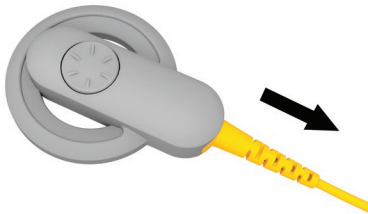

- 2 Use your fingernails to hold the coil cable grip and pull the coil cable firmly off the coil. Do not tug on the flexible part of the coil cable.



- 3 Push the new coil cable into the processor until it clicks into place.

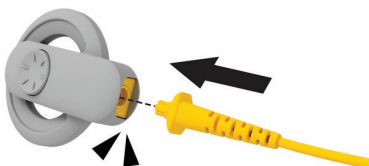

- 4 Push the new coil cable into the coil until it clicks into place.



CARE

Change the coil magnet

The strength of the Cochlear™ Nucleus® CP900 Series Coil Magnet can be adjusted by changing the coil magnet depth or by replacing with a coil magnet of a different strength. If the coil magnet is too weak or strong the coil may fall off or cause discomfort.



WARNING


If the coil magnet is too strong or is in contact with your skin, you may develop a pressure sore. If you experience any discomfort, contact your clinician.

NOTE

If your clinician gives you a Cochlear™ Nucleus® CP900 Series Coil Magnet Reverse Polarity, use it as described for a normal magnet.



Change the strength of the coil magnet

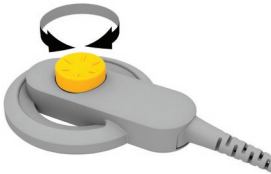

- 1 With finger grips facing up, turn the coil magnet clockwise to increase the coil magnet's strength.
- 2 Turn the coil magnet anticlockwise to reduce its strength.


NOTE

If your clinician uses the strongest magnet (a Cochlear™ Nucleus® CP900 Series Coil (6M)), it is part of the coil and cannot be adjusted like other magnets.



See *Processors and parts* on page 83.


Change the coil magnet

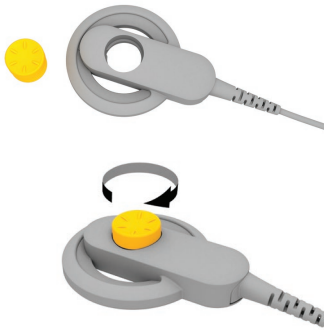

- 1 With finger grips facing up, unscrew the coil magnet anticlockwise and remove.
- 2 Insert the new coil magnet into the coil hole and turn it clockwise until it is at a comfortable depth.

Use a coil spacer

If changing to a weaker magnet does not resolve discomfort, or coil performance is affected by your coil being too close to the implant, you can use a Cochlear™ Nucleus® Coil Spacer.

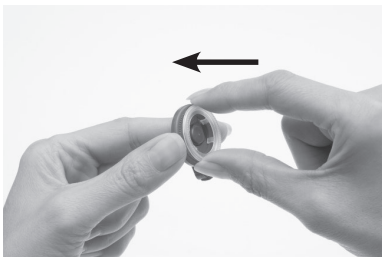

- 1 Hold the coil between the thumb and index finger with the underside of the coil facing up.
- 2 Hold the coil spacer in the other hand and align the coil spacer clips with the grooves in the coil.
- 3 Clip the coil spacer into place. Do not force the coil spacer into place. It should fit easily.



Care for your acoustic component

The acoustic component is used only by people who are set up to use their processor in Hybrid mode.

If you use Power or Plus Domes

Every day

- Check the dome and its filter. Change the dome whenever its filter is clogged with wax, or whenever the dome is dirty, discoloured or losing shape.
- When changing the dome, also check the wax filter in the speaker unit. Change the speaker unit filter whenever it is clogged with wax or if the sound from the speaker unit does not sound normal.



CAUTION

Handle and store the acoustic component very carefully. Do not drop the acoustic component. Do not twist the speaker unit cable. Keep the speaker unit away from the coil magnet.

If you use RITE Power or Power Micro Moulds

Every day

- Clean the earmould and speaker unit cable with the supplied brush or a soft dry cloth. Wipe any traces of wax away from the opening of the vent.
- Check if the vent is clogged, and use the supplied vent cleaner to remove wax as needed. Do NOT use the vent cleaner to clean any other part.
- Check the wax filter and change whenever it is clogged with wax or if the sound from the speaker unit does not sound normal.


NOTE

Power Micro Moulds can use Cochlear EAC200 Series WaxStop Wax Management System filters ONLY.


RITE Power Moulds can use either Cochlear EAC200 Series NoWax Wax Management System or WaxStop filters.




Replace Power Domes

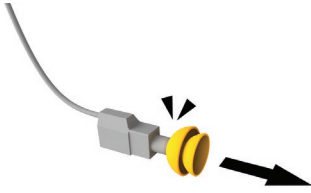

- 1 Pull the dome off using your fingernails if necessary and discard.

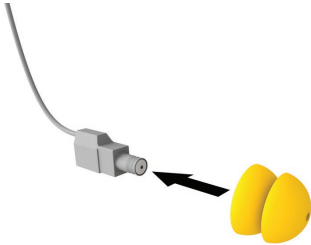

- 2 Push the new dome directly onto the speaker unit.

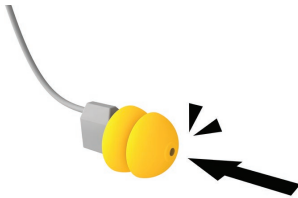

- 3 Click it securely into place.

Replace Plus Domes

Follow the directions for replacing Power Domes and add the following step:



- 4 When fitting a Plus Dome, the outer sleeve should face down.

See *Acoustic component* on page 84.

Change wax filters


Replace a NoWax filter

RITE Power Moulds (only) can use a Cochlear EAC200 Series NoWax Wax Management System filter in the sound outlet. Change the filter whenever it is clogged with wax or whenever the processor does not sound normal.




CAUTION


Wax filters must ALWAYS be used with the acoustic component and replaced regularly.

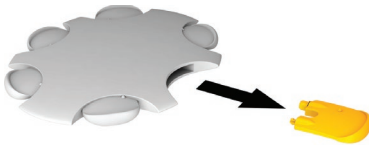

- 1 Remove a new NoWax tool from its shell. One point on the tool has a removal tool and the other point has the new filter.

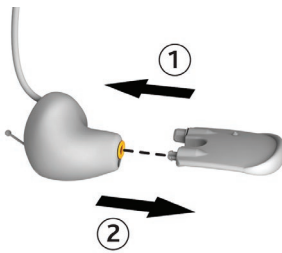

- 2 Insert the removal tool into the existing wax filter and work it out of the RITE Power Mould.

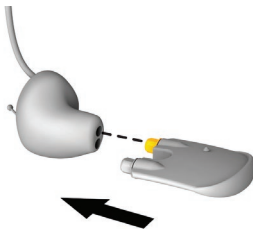

- 3 Push the new filter into the sound outlet on the earmould and dispose of the tool and used filter.




CARE

Replace a WaxStop filter


Power Micro Moulds use a Cochlear EAC200 Series WaxStop Wax Management System filter in the sound outlet. RITE Power Moulds use either a WaxStop filter or a NoWax filter. Power and Plus Domes must be worn with a WaxStop filter in the speaker unit.




- 1 Remove a WaxStop tool from its shell. One end of the tool has a new filter and the other has a removal tool.

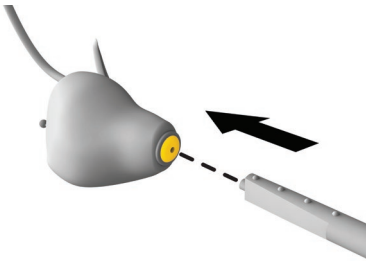

- 2 Insert the removal tool into the existing wax filter.

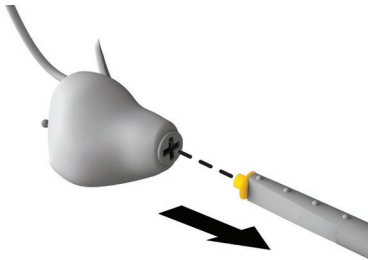

- 3 Pull the used filter out slowly, keeping the tool straight.

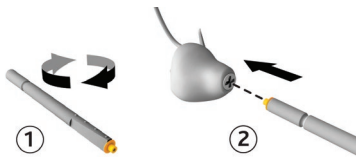

- 4 Flip the tool and push the new filter into the sound outlet. Carefully pull the tool free and dispose of the used tool and filter.



CARE








Lights and Beeps

Lights	68
Beeps	70




Lights

Your clinician can set up your processor to show some or all of the following light indications.

Turning on and off



LIGHT	WHAT IT MEANS
 ... Quick blue flashes	Processor flashes while receiving sound from telecoil/audio accessory (if set up by your clinician*).
 ... Quick green flashes	Processor flashes while receiving sound from microphones (if set up by your clinician*).
    Quick green flashes	Turning on and changing programs. Number of flashes indicates the number of the current program.
 Steady orange	Turning off processor.

Locking buttons




LIGHT	WHAT IT MEANS
 Green flash followed by orange	Locking processor buttons.
 Orange flash followed by green	Unlocking processor buttons.
 Orange flash when pressing buttons	Processor buttons are locked.






Wireless Accessories

LIGHT	WHAT IT MEANS
 Quick blue flash	Processor flashes when pairing to Wireless Accessory is successful.
 Quick blue flashes	Processor flashes while receiving sound from Wireless Accessory.

Telecoil/plug-in audio accessories

LIGHT	WHAT IT MEANS
 Quick blue flashes	Processor flashes while receiving sound from telecoil/plug-in audio accessory (if set up by your clinician*).
 Long flash of blue	Changing between microphones and telecoil/plug-in audio accessory.
 Long flash of green	Changing between the telecoil/plug-in audio accessory and microphones.

Alerts




LIGHT	WHAT IT MEANS
 Flash of orange every second	Processor flashes while coil is off (or connected to the wrong implant).
 Orange flashes	Processor battery is empty. Charge battery.
 Steady orange	Fault. Contact your clinician. Stays on until the issue is resolved.

* See your Cochlear™ Nucleus® CR230 Remote Assistant User Guide for details about modes.




Beeps

Your clinician can set up your processor so you can hear the following beeps. The beeps are only audible to the recipient.



Turning on and off

BEEP	WHAT IT MEANS
 Short high beeps	Changing the program. The number of beeps indicates the number of the selected program.
 Short high beep	Changing volume or sensitivity level.
 Short high then short low beep	When changing volume or sensitivity, indicates upper or lower limit of volume/sensitivity reached.


Locking buttons

BEEP	WHAT IT MEANS
 Short low beep when you press a button	Processor buttons are locked.
 Short high then short low beep	Locking processor buttons.
 Short low then short high beep	Unlocking processor buttons.




Wireless Accessories

BEEP	WHAT IT MEANS
 3-tone chime	Connecting with Wireless Accessory to begin streaming
 Short beep	When stopping streaming audio.




Telecoil/plug-in audio accessories

BEEP	WHAT IT MEANS
 Long high beep	Switching between using the microphones and the telecoil or plug-in audio accessory.

Alerts

BEEP	WHAT IT MEANS
 2 Short low beeps	Processor battery is low. Recharge or replace battery.
 Short low beeps for 4 seconds	Battery is empty and processor is turning off. Recharge or replace battery.
 4 long low beeps over 4 seconds	General fault. Consult your clinician.

Adjusting bass and treble (Remote Assistant only)

BEEP	WHAT IT MEANS
 Loud medium beep	Adjusting master volume level.
 Loud long high beep	Adjusting treble level.
 Loud long low beep	Adjusting bass level.




Advanced mode

Advanced mode


Sound processors can be set up in Advanced mode by your clinician to give you extra control over the volume and sensitivity of the sound you hear. Advanced mode uses a different set of button presses.



Turn processor on
Short-press the lower button.

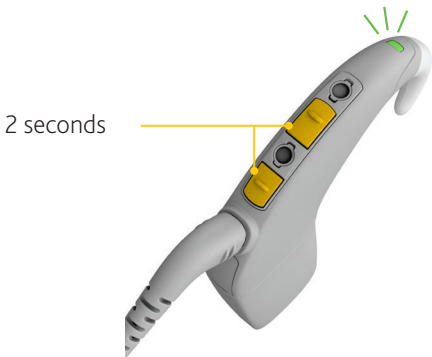


Turn processor off
Long-press upper and lower buttons at the same time.

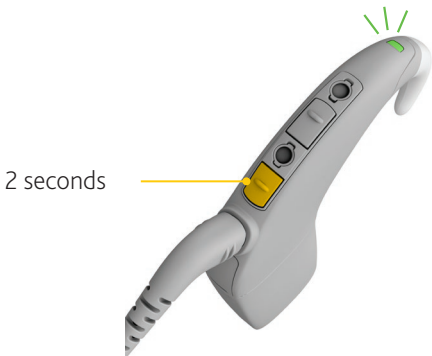

Change program
Long-press the lower button.


Turn telecoil/audio accessory on and off

Long-press the upper button to switch between microphones and telecoil/audio accessory.

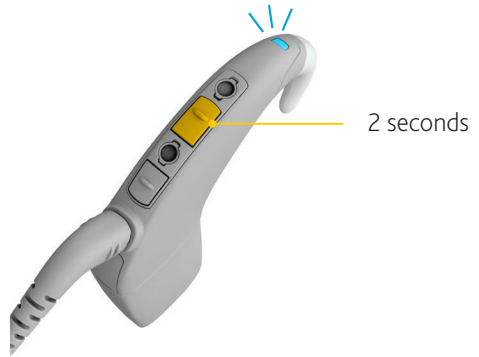

Lock/unlock processor buttons

Short-press both the upper and lower buttons at the same time.

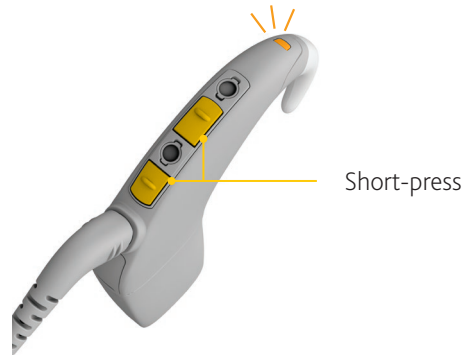


Change volume or sensitivity

Choose between volume and sensitivity (if available) and adjust them with your remote.

You can also use the buttons on your sound processor to adjust volume or sensitivity. Volume controls the level of sound you hear. Sensitivity controls the range of sounds picked up by the processor e.g. soft, background, close or far away.

You can select from up to 10 volume levels or 21 sensitivity levels.

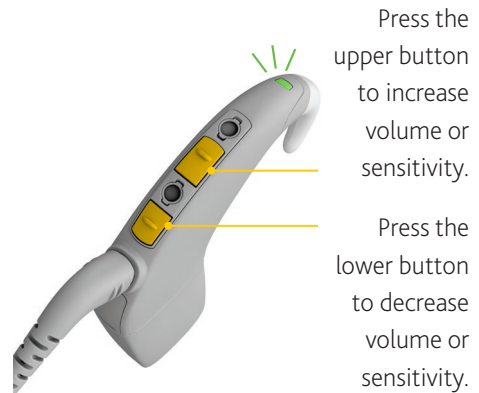

Troubleshoot

PROBLEM	RESOLUTION
You experience tightness, discomfort or develop a skin irritation at your implant site.	<ol style="list-style-type: none">1. Your coil magnet may be too strong or in contact with your skin. Change to a weaker magnet. See <i>Change the coil magnet</i> on page 58.2. If the problem continues, contact your clinician.
The process or coil become hot.	<ol style="list-style-type: none">1. Remove the processor and coil from your head immediately and contact your clinician.
You are not sure what the processor beeps or indicator light flashes mean (on top of the processor).	See <i>Lights and Beeps</i> on page 67.
You want to perform a regular check on your processor	See <i>Care for your sound processor</i> on page 46.
You want to confirm your processor is receiving sound.	<ol style="list-style-type: none">1. Check the light on the top of the processor. See <i>Lights and Beeps</i> on page 68.2. If you use a CR230 Remote Assistant, check the sound meter on the status screen.3. If you have monitor earphones, a hearing person can listen to the sound received by the processor. See <i>Plug-in audio accessories</i> on page 20 or the Accessories Supplement.

PROBLEM	RESOLUTION
<p>Processor will not turn on/ buttons will not respond.</p>	<ol style="list-style-type: none"> 1. Try turning the processor on again. See <i>Turn on and off</i> on page 14. 2. Try unlocking the buttons. See <i>Lock and unlock buttons</i> on page 18. 3. Recharge/replace the battery. See <i>Replace the battery</i> on page 7 and <i>Charge rechargeable battery modules</i> on page 10. 4. Check the battery contacts are free of dirt and dust. See <i>Water, sand and dirt</i> on page 48. 5. If the problem continues, contact your clinician.
<p>You do not hear sound or sound is intermittent.</p>	<ol style="list-style-type: none"> 1. If you use Advanced mode, turn up volume or sensitivity. See <i>Change volume or sensitivity</i> on page 75. 2. Try a different program. See <i>Change programs</i> on page 15. 3. Recharge/replace the battery. See <i>Replace the battery</i> on page 7 and <i>Charge rechargeable battery modules</i> on page 10. 4. Make sure the coil cable is fully inserted into the socket on the coil. 5. If the problem continues, contact your clinician.

PROBLEM	RESOLUTION
<p>Sound is too quiet or muffled.</p>	<ol style="list-style-type: none"> 1. Try changing the microphone protectors. See <i>Replace microphone protectors</i> on page 49. 2. If you use Advanced mode, turn up volume or sensitivity. See <i>Change volume or sensitivity</i> on page 75. 3. Try a different program. See <i>Change programs</i> on page 15. 4. If you use a CR230 Remote Assistant, see the user guide for further troubleshooting. 5. If the problem continues, contact your clinician.
<p>Sound is too loud or uncomfortable.</p>	<ol style="list-style-type: none"> 1. If you use Advanced mode, turn down volume or sensitivity. See <i>Change volume or sensitivity</i> on page 75. 2. Try a different program. See <i>Change programs</i> on page 15. 3. If the problem continues, contact your clinician.
<p>You hear intermittent sound, a buzzing sound or distorted speech.</p>	<ol style="list-style-type: none"> 1. Check for sources of interference such as radio and TV transmission towers, shopping centres, airport security systems and mobile phones. 2. Try moving away from any electronic device that may be causing interference. 3. Auto telecoil may be activated. Try turning it off. See <i>Telecoil</i> on page 23. 4. If the problem continues, contact your clinician.

PROBLEM	RESOLUTION
You do not hear sound from a Wireless Accessory.	<ol style="list-style-type: none">1. Check that the Wireless Accessory is charged and turned on.2. Check that the Wireless Accessory is paired with your processor.3. Check the volume of the Wireless Accessory.4. If you use a CR230 Remote Assistant, use the sound meter to check the processor is receiving sound from the accessory.5. If you use a CR230 Remote Assistant, check and adjust the accessory/microphone mixing ratio.6. If available, try a different processor.
You do not hear sound from a plug-in audio accessory.	<ol style="list-style-type: none">1. Check that the accessory cable is fully inserted into the socket on the processing unit.2. Check that the other end of the accessory cable is fully inserted into the audio source.3. Check that the audio source is turned on and working properly.4. Check the volume on the audio source.5. If you use a CR230 Remote Assistant, use the sound meter to check the processor is receiving sound from the accessory.6. If you use a CR230 Remote Assistant, check and adjust the accessory/microphone mixing ratio.7. If available, try a different processor.

PROBLEM	RESOLUTION
Your processor does not automatically turn on the telecoil.	<ol style="list-style-type: none"> 1. If you use a CR210 Remote Control or CR230 Remote Assistant, check that Auto telecoil is turned on. 2. Wait until you hear speech. If Auto Telecoil is turned on, it should switch from standby to receiving when it detects speech. 3. Check with your clinician if Auto Telecoil is enabled on your processor. 4. Manually turn on the telecoil by pressing the upper button on your processor. See <i>Telecoil</i> on page 23.
You have accidentally changed processor settings.	<ol style="list-style-type: none"> 1. If you use a CR230 Remote Assistant, reset your processor's settings. 2. If you do not use a CR230 Remote Assistant, contact your clinician.
The processor gets wet.	<ol style="list-style-type: none"> 1. Dry the processor with a soft cloth, change the microphone protectors and place the processor in the dry aid kit provided by Cochlear for 8 hours. See <i>Water, sand and dirt</i> on page 48.
The processor switches off automatically.	<ol style="list-style-type: none"> 1. Ask your clinician if "Auto power off" is enabled. If it is, then the processor will switch off when not connected to the implant for more than two minutes. 2. Change the battery/batteries. See <i>Replace the battery</i> on page 7 or <i>Change disposable batteries</i> on page 11.

PROBLEM	RESOLUTION
The processor will not turn off.	<ol style="list-style-type: none"> 1. Check you are pressing the correct button. See <i>Turn on and off</i> on page 14. 2. Check the processor is not locked. See <i>Lock and unlock buttons</i> on page 18. 3. Detach the battery module from the processing unit. See <i>Remove the battery</i> on page 7.
When using the LiteWear cable, the processor does not work.	<ol style="list-style-type: none"> 1. Check the LiteWear cable is properly attached to both the battery module and the processor. 2. Try turning on your processor again. 3. Check the battery contacts and ensure they are free of dirt, dust or corrosion.
The indicator light on the battery charger shows a fully charged battery is still charging.	<p>This won't cause any damage to the battery as the charge cycle on a fully charged battery is very short.</p>
The indicator light on the battery charger is flashing orange	<ol style="list-style-type: none"> 1. Check that only rechargeable battery modules are connected to the battery charger. 2. Replace rechargeable battery module.
Disposable batteries are not lasting as long as usual.	<ol style="list-style-type: none"> 1. Clean all connections and parts of the tamper resistant battery module. 2. Clean all connections on the processing unit. 3. Try replacing the coil with a new coil. 4. If the problem continues, contact your clinician.

Warnings

NOTE

Please also refer to your *Important Information* document for essential advice that applies to Cochlear implant systems.



For parents and carers

- Magnetic resonance imaging (MRI) is contraindicated except under the circumstances described in your *Important Information* document.
- Remove the sound processor before entering a room where an MRI scanner is located.
- Sound processor covers and coil covers may not be suitable for children. Removable parts can be lost or may be a choking hazard.
- Keep chemicals from dry aid kits away from small children. Swallowing this material can cause serious internal injuries.
- Removable parts of the system (e.g. Snugfit, Mic Lock) can be lost or may be a choking hazard.
- Parents and carers are advised that unsupervised use of long cables (e.g. coil or accessory cables) may present a risk of strangulation.
- Make sure all cables used by a child are securely attached to their clothing.
- Remove the processor and coil immediately if there is any discomfort or pain (e.g. if device becomes hot, or sound is uncomfortably loud) and inform clinician. Carer must routinely check the device for signs of overheating and for signs of discomfort or pressure sores at the coil site.
- Do not allow children to replace batteries without adult supervision.
- Do not touch the battery charger contacts or allow children to use the battery charger without adult supervision.

- Keep all parts of the Cochlear Nucleus CP900 Series Tamper Resistant Battery Module and button cell batteries out of reach of children. Removable parts can be lost or may be a choking hazard.
- Contact your clinician for advice on choosing the most suitable acoustic component for children.
- Do not allow children to replace an acoustic component without adult supervision.
- Keep all parts of the acoustic component out of reach of children. Removable parts (e.g. Power/Plus domes and wax management systems) can be lost or may be a choking hazard.
- Carers should routinely check that the acoustic component is working at a comfortable volume level. See your clinician if the sound level ever becomes uncomfortable.
- The tamper resistant earhook fitting tool and fitting pin are sharp objects. They are not recommended for use by children.

Processors and parts

- No modification of this equipment is allowed.
- Remove the monitor earphones immediately if the sound level is uncomfortably loud and inform your clinician.
- If the coil magnet is too strong or is in contact with the skin, pressure sores may develop at the coil site. If this happens, or if you experience any discomfort in this area, contact your clinician.
- Do not push the volume too high for comfort in case a loud noise occurs nearby.
- Do not place the device or accessories inside any part of your body (e.g. nose, mouth).
- Do not let your processor or parts entangle with any jewellery (e.g. earhook and earrings) or machinery.

- Do not apply continued pressure to the coil when in contact with the skin (e.g. sleeping while lying on coil, or using tight fitting headwear). Do not attach or wear a body-worn battery configuration beneath layers of clothing in direct contact with the skin (e.g. in pocket or under a child's vest).
- Do not place the processor or parts in any household devices (e.g. microwave oven, dryer).
- Do not modify your processor. Warranty will be void if modified.
- Do not use a dry aid kit that has an Ultra Violet C (UVC) lamp (e.g. do not use the Freedom Dry and Store).
- Your device contains magnets that should be kept away from life supporting devices (e.g. cardiac pacemakers and ICDs (implantable cardioverter defibrillators) and magnetic ventricular shunts), as the magnets may affect the function of these devices. Keep your processor at least 15 cm (6 in) from such devices. Contact the manufacturer of the specific device to find out more.
- Store spare coil magnets safely and away from cards that may have a magnetic strip (e.g. credit cards, bus tickets, etc).
- Advanced mode users – if you find you need to adjust the volume often, or if adjusting volume ever causes discomfort, consult your clinician.

Acoustic component

- Refit and replace domes regularly to ensure optimal fitting and retention. Consult your clinician for advice on acoustic component misfits or signs of discomfort.
- Remove acoustic component if there is any discomfort or pain (e.g. sound level is uncomfortably loud) and inform your clinician.
- If the Power Dome or Plus Dome is not securely attached to your acoustic component, it could fall off in your ear canal. If this happens, see your medical practitioner.

Batteries

- Dispose of used batteries promptly and carefully, in accordance with local regulations. Keep away from children.
- Do not recharge disposable batteries.
- Do not disassemble, deform, immerse in water or dispose of batteries in fire.
- Replace the whole set of two disposable batteries with high power Cochlear supplied or recommended zinc air batteries only. We recommend 675 (PR44) zinc air batteries designed for cochlear implant use. We do not recommend using silver oxide or alkaline batteries.
- Do not mix disposable batteries that differ by manufacturer, brand, type, age or previous usage.
- Insert batteries correctly.
- Only use rechargeable batteries and battery chargers supplied or recommended by Cochlear. Use of other batteries or battery chargers may result in harm or injury.
- Charge rechargeable batteries before use.
- Do not short-circuit batteries (e.g. do not let terminals of batteries contact each other, do not place batteries loose in pockets, etc.).
- Store unused batteries in original packaging, in a clean and dry place. When processor is not in use, remove the disposable or rechargeable batteries and store separately in a clean and dry place. Wipe batteries with a clean dry cloth if they become dirty.
- Do not expose batteries to heat (e.g. never leave batteries in sunlight, behind a window or in a car).
- Do not use damaged or deformed batteries. If skin or eyes come into contact with battery fluid or liquid, wash out with water and seek medical attention immediately.
- Never put batteries in mouth. If swallowed, contact your physician or local poison information service.

Other information

Physical configuration

The processing unit comprises:

- Two omni-directional microphones for receiving sound.
- An internal telecoil for receiving magnetic fields radiated by phones, neckloops and roomloops.
- Custom analogue and digital integrated circuits with digital signal processing (DSP) and bi-directional wireless communication capabilities.
- A tri-colour visual indication of processor function or problem.
- Two buttons allow user control of key features.
- CP910 only – custom 4-pin accessory connector for connection
- Custom 4-pin coil connector for connection of the coil cable.
- A range of earhooks.

The batteries provide power to the processor. The coil acts as a transformer coupling that transfers energy and data to the implant.

Wireless communication link

The remote control/remote assistant wireless communication link operates in the 2.4 GHz ISM band using GFSK (Gaussian frequency shift keying) on 5 channels. The link uses a proprietary bi-directional communication protocol and operates over a distance of up to 2 metres from the processor. When interference is present, the wireless communication link switches between the 5 channels to find a channel where the interference least affects the operation of the link. The remotes indicate via their displays when the processor is not within operating distance, and when the link has been interrupted due to interference (see the relevant remote user guide for more information).

Materials

- Processing unit: copolyester
- Battery modules (all types) are made of copolyester.
- Coil is made of polypropylene (PP), thermoplastic elastomer (TPE).
- Coil magnet casing is made of acrylonitrile butadiene styrene (ABS).
- Coil cable sheath is made of polyvinyl chloride (PVC).
- Coil cable plugs are made of PP & TPE.
- LiteWear cable is made of polyester, TPE & PVC (same construction as the coil cable) sheath is PVC; plugs/strain relief are PP and TPE.
- LiteWear connection shoes are made of copolyester.
- LiteWear retention case is made of polyamide.
- LiteWear Fixing Aids are made of polycarbonate.
- Processor covers are made of polystyrene and polycarbonate.
- Coil covers are made of polycarbonate.
- Hybrid materials for outer sheathing of the housing adaptor and the outer cable sheathing is made of polyamide.
- Custom earmoulds are made of Stereo Lithography Acrylic covered with a UV-cured lacquer.
- Plus Domes are made of elastosil silicone polymer.
- Power Domes are made of TPE.

Battery life, charge cycles and lifespan

- Battery life means the time a device will run before the disposable batteries must be replaced, or the rechargeable batteries recharged.
- Battery charge cycle is a full charge and discharge of the rechargeable battery.
- Battery lifespan means the total number of charge cycles a rechargeable battery will last before the battery life degrades to 80% of its original fully-charged capacity.

Product component dimensions *(Typical values)*

COMPONENT	LENGTH	WIDTH	DEPTH	DIAMETER
CP910 processing unit with medium earhook and standard rechargeable battery module	51.3 mm	9.0 mm	47.7 mm	N/A
CP920 processing unit with medium earhook and compact rechargeable battery module	39.1 mm	9.0 mm	47.7 mm	N/A
Standard tamper resistant battery module	28.0 mm	9.0 mm	19.0 mm	N/A
Compact rechargeable battery module	19.0 mm	9.0 mm	19.0 mm	N/A
Standard rechargeable battery module	28.0 mm	9.0 mm	19.0 mm	N/A
Coil	N/A	N/A	8.0 mm	30.8 mm

Weight *(Typical values. All weights are measured with a medium earhook.)*

COMPONENT	WEIGHT
CP910 processing unit (no battery module)	5.7g
CP920 processing unit (no battery module)	5.0g
CP910 processing unit with compact rechargeable battery module	10.5g
CP910 processing unit with standard rechargeable battery module	13.0g
CP910 processing unit with standard tamper resistant battery module (including two 675 zinc air batteries)	12.8g
CP920 processing unit with compact rechargeable battery module	9.8g
CP920 processing unit with standard rechargeable battery module	12.4g
CP920 processing unit with standard tamper resistant battery module (including two 675 zinc air batteries)	12.1g
Coil (without coil magnet)	4.5g

Operating characteristics

Processing unit

CHARACTERISTIC	VALUE/RANGE
Sound input frequency range	100 Hz to 8 kHz
Wireless technology	Proprietary low power bi-directional wireless link
RF frequency	2.4 GHz
Operating voltage	2.00 V to 4.25 V
Power consumption	20 mW to 100 mW
Charge cycles	≥ 80% capacity after 365 charge / discharge cycles at room temperature
Button functions	Turn processor on and off, turn telecoil/audio accessory on and off, change program, lock/unlock buttons, change sensitivity or volume level.
Wireless transmission range	Up to 2 m

Battery module

TYPE	CAPACITY/VOLTAGE RANGE
Standard tamper resistant battery module	Two PR44 (zinc air) button cell batteries. 1.45 V (nominal) each. Cochlear recommends 675 zinc air batteries designed for cochlear implant use.
Compact rechargeable battery module	120 mAh/ 3.0 V to 4.2 V
Standard rechargeable battery module	205 mAh/ 3.0 V to 4.2 V

Coil

CHARACTERISTIC	VALUE/RANGE
Operating voltage	2.0 V to 2.6 V
Operating frequency (standard coil)	5 MHz
Operating frequency (N22 coil)	2.5 MHz

Environmental conditions

CONDITION	MINIMUM	MAXIMUM
Storage & transport temperature	-10°C (14°F)	+55°C (131°F)
Storage & transport humidity	0% RH	90% RH
Operating temperature	+5°C (41°F)	+40°C (104°F)
Operating relative humidity	0% RH	90% RH
Operating pressure	700 hPa	1060 hPa

Environmental protection

Your sound processor contains electronic components subject to the Directive 2002/96/EC on waste electrical and electronic equipment.

Help protect the environment by not disposing of your sound processor or batteries with your unsorted household waste. Please recycle your sound processor according to your local regulations.

Standard tamper resistant battery module

Check the battery manufacturer's recommended operating conditions for disposable batteries used in your processor.

Certification and applied standards

Your device fulfils the essential requirements listed in Annex 1 of the EC directive 90/385/EEC on Active Implantable Medical Devices as per the conformity assessment procedure in Annex 2.



The year in which authorisation to affix the CE mark was granted was 2015.

It also fulfils the essential requirements of Directive 2014/53/EU on Radio Equipment as per the conformity assessment procedure in Annex II. The declaration of conformity may be consulted at www.cochlear.com/wps/wcm/connect/intl/about/company-information/declaration-of-conformity



Equipment classification

Your sound processor is internally powered equipment Type B applied part as described in the international standard IEC 60601-1:2005/A1:2012, Medical Electrical Equipment – Part 1: General Requirements for Basic Safety and Essential Performance.

FCC (Federal Communications Commission) and Canadian IC compliance

This device complies with part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- This device may not cause harmful interference.
- This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications made to this equipment not expressly approved by Cochlear Limited may void the FCC authorisation to operate this equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can

be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet or a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

FCC ID: WTOP900

IC ID: 8039A-P900

Radio compliance for Japan

This device is granted pursuant to the Japanese Radio Law (電波法) and the Japanese Telecommunications Business Law (電気通信事業法).

This device should not be modified (otherwise the granted designation number will become invalid).



Radio compliance for Korea

Warning: This radio equipment has the possibility of radio interference during operation.



Labelling symbols

The following symbols may appear on your processor or remote components and/or packaging:



Refer to instruction manual



Specific warnings or precautions associated with the device, which are not otherwise found on the label



Manufacturer



Authorised representative in the European Community



Catalogue number



Serial number



Batch code



Date of manufacture



Temperature limits



CE registration mark with notified body number



Radio compliance certification for Australia and New Zealand



Radio compliance certification for Japan



Radio compliance certification for Korea



Do not re-use

Rx Only

By prescription



Recyclable material



Coil guide. Used to position the coil against the remote control when pairing.



Dispose of electrical components in accordance with your local regulations



Dispose of lithium-ion batteries in accordance with your local regulations



Type B applied part

Ingress Protection Rating

- Protected against access of solid foreign objects greater than or equal to 1.0mm diameter.
- Protected against failure from dust penetration.
- Protected against failure from splashing water.

IP54

Ingress Protection Rating

- Protected against access of solid foreign objects greater than or equal to 1.0mm diameter.
- Protected against failure from dust penetration.
- Protected against failure from temporary immersion in water.

IP57

Legal statement

The statements made in this guide are believed to be true and correct as of the date of publication. However, specifications are subject to change without notice.

© Cochlear Limited 2017

Cochlear™



Cochlear Ltd (ABN 96 002 618 073) 1 University Avenue, Macquarie University, NSW 2109, Australia
Tel: +61 2 9428 6555 Fax: +61 2 9428 6352

Cochlear Ltd (ABN 96 002 618 073) 14 Mars Road, Lane Cove, NSW 2066, Australia
Tel: +61 2 9428 6555 Fax: +61 2 9428 6352

COCHLEAR Cochlear Deutschland GmbH & Co. KG Karl-Wiechert-Allee 76A, 30625 Hannover, Germany
Tel: +49 511 542 770 Fax: +49 511 542 7770

Cochlear Americas 13059 E Peakview Avenue, Centennial, CO 80111, USA
Tel: +1 303 790 9010 Fax: +1 303 792 9025

Cochlear Canada Inc 2500-120 Adelaide Street West, Toronto, ON M5H 1T1, Canada
Tel: +1 416 972 5082 Fax: +1 416 972 5083

Cochlear AG EMEA Headquarters, Peter Merian-Weg 4, 4052 Basel, Switzerland
Tel: +41 61 205 8204 Fax: +41 61 205 8205

Cochlear Europe Ltd 6 Dashwood Lang Road, Bourne Business Park, Addlestone, Surrey KT15 2HJ, United Kingdom
Tel: +44 1932 26 3400 Fax: +44 1932 26 3426

Cochlear Benelux NV Schaliënhoedreef 20 i, B-2800 Mechelen, Belgium
Tel: +32 15 79 55 11 Fax: +32 15 79 55 70

Cochlear France S.A.S. 135 Route de Saint-Simon, 31035 Toulouse, France
Tel: +33 5 34 63 85 85 (International) or 0805 200 016 (National) Fax: +33 5 34 63 85 80

Cochlear Italia S.r.l. Via Larga 33, 40138 Bologna, Italy
Tel: +39 051 601 53 11 Fax: +39 051 39 20 62

Cochlear Nordic AB Konstruktionsvägen 14, 435 33 Mölnlycke, Sweden
Tel: +46 31 335 14 61 Fax: +46 31 335 14 60

Cochlear Tibbi Cihazlar ve Sağlık Hizmetleri Ltd. Şti.

Çubuklu Mah. Boğaziçi Cad., Boğaziçi Plaza No: 6/1, Kavacık, TR-34805 Beykoz-Istanbul, Turkey
Tel: +90 216 538 5900 Fax: +90 216 538 5919

Cochlear (HK) Limited Room 1204, 12/F, CRE Building, No 303 Hennessy Road, Wanchai, Hong Kong SAR
Tel: +852 2530 5773 Fax: +852 2530 5183

Cochlear Korea Ltd 1st floor, Cheongwon Building 33, Teheran-ro 8 gil, Gangnam-gu, Seoul, Korea
Tel: +82 2 533 4450 Fax: +82 2 533 8408

Cochlear Medical Device (Beijing) Co., Ltd

Unit 2208-2212, Tower B, Gerndale Building, 91 Jianguo Road, Chaoyang District, Beijing 100022, P.R. China
Tel: +86 10 5909 7800 Fax: +86 10 5909 7900

Cochlear Medical Device Company India Pvt. Ltd.

Ground Floor, Platina Building, Plot No C-59, G-Block, Bandra Kurla Complex, Bandra (E), Mumbai – 400 051, India
Tel: +91 22 6112 1111 Fax: +91 22 6112 1100

株式会社日本コクレア (Nihon Cochlear Co Ltd) 〒113-0033 東京都文京区本郷2-3-7 お茶の水元町ビル
Tel: +81 3 3817 0241 Fax: +81 3 3817 0245

Cochlear Middle East FZ-LLC

Dubai Healthcare City, Al Razi Building 64, Block A, Ground Floor, Offices IR1 and IR2, Dubai, United Arab Emirates
Tel: +971 4 818 4400 Fax: +971 4 361 8925

Cochlear Latinoamérica S.A.

International Business Park, Building 3835, Office 403, Panama Pacifico, Panama
Tel: +507 830 6220 Fax: +507 830 6218

Cochlear NZ Limited

Level 4, Takapuna Towers, 19-21 Como St, Takapuna, Auckland 0622, New Zealand
Tel: +64 9 914 1983 Fax: 0800 886 036

www.cochlear.com

ACE, Advance Off-Stylet, AOS, AutoNRT, Autosensitivity, Beam, Button, Carina, Cochlear, Cochlear SoftWear, コクレア,
Codacs, Contour, Contour Advance, Custom Sound, ESprit, Freedom, Hear now. And always, Hugfit, Hybrid, inHear, Invisible
Hearing, Kanso, MET, MicroDrive, MP3000, myCochlear, mySmartSound, NRT, Nucleus, 科利耳, Off-Stylet, SmartSound,
Softip, SPrint, True Wireless, the elliptical logo and Whisper are either trademarks or registered trademarks of Cochlear Limited.
Ardium, Baha, Baha SoftWear, BCDrive, Dermalock, EveryWear, Vistafix and WindShield are either trademarks or registered
trademarks of Cochlear Bone Anchored Solutions AB. Bluetooth is a registered trademark of Bluetooth SIG.
© Cochlear Limited 2017

Hear now. And always

390443 ISS12 OCT17