



Cochlear[™] Osia[®] OSI200 Implant

Magnetic Resonance Imaging (MRI) Guidelines

About this guide

This guide applies to Cochlear™ Osia® OSI200 Implant. It is intended for:

- Specialised health care professionals who prepare and perform MR scans
- Physicians who refer a Cochlear Osia implant recipient for an MR scan
- Cochlear Osia implant recipients and/or their carers.

This guide provides information about the safe application of an MR scan on Cochlear Osia OSI200 Implant recipients.

MR scans performed under different conditions than those presented in this guide may result in severe patient injury or device malfunction.

Due to the risks associated with using MRI with an implanted medical device, it is important to read, understand, and comply with these instructions to prevent potential harm to the patient and/or device malfunction.

This guide should be read in conjunction with the relevant documents that accompany a Cochlear Osia OSI200 Implant, such as the *Physician's Guide and Important information for Osia implant recipients*.

For more information, contact Cochlear by calling your regional Cochlear office – contact numbers are available on the back cover of this guide or visit www.cochlear.com/mri

If you are a consumer, please seek advice from your medical practitioner or health professional prior to an MR scan.

Symbols used in this guide



NOTE

Important information or advice.



CAUTION (NO HARM)

Special care to be taken to ensure safety and effectiveness.
Could cause damage to equipment.



WARNING (HARMFUL)

Potential safety hazards and serious adverse reactions.
Could cause harm to person.

Contents

About this guide	2	Identifying the Cochlear Osia implant 10	
Symbols used in this guide	2	X-ray information	10
Preparation prior to an MRI examination	4	X-ray guidelines	11
Cooperation between specialists	4	Performing MRI safely	12
Determine eligibility for MRI	5	Implant magnet conditions for MRI	12
Risks associated with MRI and Cochlear Osia implants	6	Scan conditions and SAR limits	12
Considerations for implant magnet removal	7	Image interference and artefacts	16
Preparation for conducting the MRI examination	8	Considerations after an MRI examination	24
Confirm the following prior to scanning:	8	With the implant magnet in place	24
Bilateral recipients	8	With the implant magnet removed	24
Performing an MR scan on other body locations	9	Disposal	24
Patient positioning	9		
Patient comfort	9		

Preparation prior to an MRI examination

These guidelines are specific to Cochlear Osia implants and supplement other MRI examination considerations specified by the MRI machine manufacturer or protocols at the MRI facility.

WARNING

Non-clinical testing has demonstrated that the Cochlear Osia OSI200 Implant, in combination with the BI300 Implant, are MRI Conditional. Patients can be scanned at 1.5 T with the magnet in place or removed. Patients can be scanned at 3 T only if the magnet has been removed. A patient with these devices can be safely scanned under the conditions described in section *“Performing MRI safely”* on **page 12 to 15**. Failure to follow these conditions may result in injury to the patient.

Cooperation between specialists

Preparing for and conducting an MRI examination for implant recipients requires cooperation between a specialist for the device and/or Osia implant physician, referring physician and radiologist/MR technologist.

Cochlear Osia implant device specialist

Knows the implant type and where to find the correct MR parameters for the implant.

Referring physician

Knows the location of the MR scan and diagnostic information required, and makes a decision on whether the implant magnet needs to be removed for the MRI examination.

Cochlear Osia implant physician

If requested by the referring physician, surgically removes the implant magnet and replaces with a non-magnetic plug. After the MR scan, the implant physician replaces it with a new sterile replacement implant magnet.

Radiologist/MR technologist

Sets up the MR scan using the correct MR parameters and counsels the implant recipient during the MRI examination.

Determine eligibility for MRI

If you are a physician referring a Cochlear Osia implant recipient for an MR scan, it is essential that you consider the following:

- Understand and inform the patient of the risks associated with MRI. See ***“Risks associated with MRI and Cochlear Osia implants”*** on **page 6**.
Also consider:
 - Timing of the implant surgery and MRI exposure.
 - Age and general health of the implant recipient and time to recover from the implant magnet surgery or potential trauma.
 - Existing or potential for tissue scarring in the location of the implant magnet.
- Understand the conditions for an MR scan and ensure that there is a clear indication for the MRI examination. See ***“Performing MRI safely”*** on **page 12**.
- The Cochlear Osia implant will create shadowing on the MR image in the vicinity of the implant, resulting in a loss of diagnostic information. Refer to the relevant artefact dimension tables in ***“Image interference and artefacts”*** on **page 16**.
- Identify if the patient has any other medical device implants, active or abandoned. If another implant is present, verify MRI compatibility before conducting an MRI examination. If MRI safety information for the implanted devices are not followed potential risks include: movement or damage to the device, weakening of the implant magnet and uncomfortable sensation or skin/tissue trauma for the patient.
- Cochlear has evaluated the interaction of implants described in this guide with other nearby implanted devices during MRI scanning, and there is no increased heating risk to the Cochlear Osia implants.
- For MR scans at 1.5 T or 3 T, identify if the implant magnet needs to be removed. See ***“Implant magnet conditions for MRI”*** on **page 12**.
- If the implant magnet needs to be removed, refer the patient to an appropriate physician to arrange for the magnet to be removed before the MR scan.
- If the implant magnet is retained for an MR scan at 1.5 T, a Cochlear Osia MRI Kit must be obtained beforehand for use during the MR scan. Contact the nearest Cochlear office or official distributor to order an MRI Kit.

Risks associated with MRI and Cochlear Osia implants

The potential risks of performing MRI examinations on patients with Cochlear Osia implants include:

Device movement

Scanning outside of the parameters contained in these guidelines may lead to the implant magnet or device moving out of position during an MRI examination causing skin/tissue trauma.

Damage to the device

MRI exposure beyond the values contained in these guidelines may cause damage to the device.

Weakening of implant magnet

- Scanning at static magnetic field strengths at values other than those contained in these guidelines may lead to a weakening of the implant magnet.
- Incorrect patient positioning prior to the MR scan or head movement during the scan may result in implant magnet demagnetisation.

Uncomfortable sensation

MRI exposure beyond the values contained in these guidelines may result in the patient perceiving sound or noise and/or pain.

Implant heating

Use the recommended SAR values contained in these guidelines to ensure the implant does not heat beyond safe levels.

Image artefact

The Cochlear Osia Implant will create shadowing on the MR image in the vicinity of the implant, resulting in a loss of diagnostic information.

If inspecting near the implant, removal of the implant magnet should be considered as MR image quality may be compromised with it in place.

Considerations for implant magnet removal

If the implant magnet needs to be removed prior to an MRI examination, close coordination is required between the specialists to perform the implant magnet removal, MR scan, and subsequent implant magnet replacement. For details on implant magnet removal, please refer to the OSI200 Implant Physician's Guide supplied with the system.

WARNING

To prevent infection, do not leave the magnet pocket empty. When removing the magnet, replace the magnet with a non-magnetic plug.

For recipients requiring multiple MRI examinations over a period of time, the implant magnet is removed and replaced with a sterile non-magnetic plug. See

"Implant magnet conditions for MRI" on page 12.

In the magnet's absence, the non-magnetic plug prevents fibrous tissue growing into the implant recess. Such growth would make implant magnet replacement difficult.

With the non-magnetic plug in place, MR scans can be safely done without the need for bandaging or use of the Cochlear Osia MRI Kit.

NOTE

While the implant magnet is removed, the recipient must wear a retainer disc to hold their sound processor coil in place. Retainer discs are available from Cochlear.

When there is no further need for MRI examinations, the non-magnetic plug is removed and replaced by a new sterile replacement implant magnet.

The sterile non-magnetic plug and sterile replacement magnet are supplied separately in sterile packs. Both are single-use items.

Preparation for conducting the MRI examination



All external components of the Cochlear Osia System (e.g. sound processors and related accessories) are MR Unsafe. The patient must remove all external components of their system before entering a room where an MRI scanner is located.

Confirm the following prior to scanning:

- The implant model has been identified. See *“Identifying the Cochlear Osia implant” on page 10.*
- For additional information for bilateral recipients, see *“Bilateral recipients” on page 8.*
- For MR scans on a body location away from the implant site, MRI safety information for the recipient’s implant model must be followed. See *“Performing an MR scan on other body locations” on page 9.*
- The implant magnet has been surgically removed when the referring physician has prescribed that the MR scan be performed with the implant magnet removed.
- If the implant magnet is retained for an MR scan at 1.5 T, a Cochlear Osia MRI Kit is available and ready for use. Refer to the Cochlear™ Osia® MRI Kit User Guide.
- Remove the sound processor before entering the MRI scan room. The sound processor is MR Unsafe.
- Position the patient to minimise discomfort. See *“Patient positioning” on page 9.*
- Discuss the sensations the recipient may experience during the MR Scan. See *“Patient comfort” on page 9.*
- Comply with the *“Scan conditions and SAR limits” on page 12.*

Bilateral recipients

CAUTION

If one of the implants is a CI22M cochlear implant without a removable magnet, MRI is contraindicated.

If a bilateral recipient has implant models (other than the CI22M cochlear implant without a removable magnet), read the MRI safety information for each implant model relevant to the recipient. Use the MRI safety information of the recipient’s implant model with the most restrictive MRI exposure requirements.

Performing an MR scan on other body locations

When an implant recipient requires an MRI on a location of their body away from the implant site, you must still follow the MRI safety information for the recipient's implant model. See *"Identifying the Cochlear Osia implant"* on **page 10** and related *"Performing MRI safely"* on **page 12**.

Patient positioning

The patient should be positioned prior to entering the MRI machine. The patient should be placed in the supine position (lying flat on back, face upward), with their head aligned with the bore axis of the MRI machine.

The patient should be advised to lie as still as possible and to not move their head during the MR scan.

CAUTION

When scanning with the implant magnet in place, ensure that the patient does not move more than 15 degrees (15°) from the centreline (Z-axis) of the bore during the MR scan.

Failure to position the patient correctly prior to the MR scan may result in increased torque on the implant and cause pain.

Patient comfort

Explain to the patient that they may sense the implant magnet moving. The MRI kit will reduce the likelihood of the implant magnet moving. However they may still sense resistance to movement as pressure on the skin. The sensation will be similar to pressing down firmly on the skin with the thumb.

If the patient experiences pain, consult the patient's physician to determine if the implant magnet should be removed or if a local anaesthetic may be applied to reduce discomfort.

CAUTION

If administering local anaesthetic, take care not to perforate the implant silicone. In addition, explain to the patient that they may perceive sounds during the MR scan. The sound levels are not hazardous.

Identifying the Cochlear Osia implant

The implant model can be found on the patient's Cochlear patient identification card. If the patient does not have their patient identification card with them, the implant type and model can be identified without surgical intervention using X-ray or the Cochlear Osia fitting software.

X-ray information

Cochlear Osia OSI200 Implants are made of metal and implanted under the skin behind the ear. Using an X-ray, the implant can be identified by its shape and the shape of the actuator unit.

Use the **Fig.1** and the **Fig.2** to assist with identifying Cochlear Osia OSI200 Implants when using an X-ray.

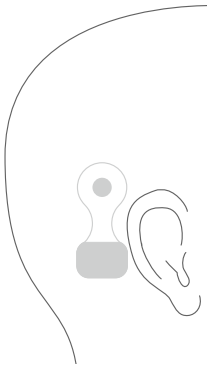


Fig.1: Approximate location of the OSI200 Implant



Fig.2: OSI200 Implant (P1170466)

X-ray guidelines

Lateral X-ray at 70 kV/ 3 mAs provides sufficient contrast to identify the implant.

A modified Stenver's view is not recommended for implant identification as implants may appear oblique.

Imaging should include an unobstructed view of antenna coils and implant bodies.

Bilateral recipients may have different implant models on either side of the head. A

lateral skull X-ray with a 15 degree cranial tube angle will offset the implants in the image, enabling identifying features to be distinguished.

Performing MRI safely

Implant magnet conditions for MRI

For some implant models and MRI field strengths, bandaging with an MRI kit is required, or the implant magnet needs to be surgically removed. Refer to the table below for information on the OSI200 Implant. Refer to the Cochlear™ Osia® MRI Kit User Guide for instructions on how to apply the MRI Kit prior to the MR scan.

Implant type	MRI field strength (T)	Remove implant magnet Yes/No	MRI Kit required Yes/No
Osia OSI200 Implant	1.5	No	Yes
	3	Yes	No

Table 1: Implant magnet conditions for MRI.

Scan conditions and SAR limits

The MRI safety information provided in these guidelines only applies to 1.5 T and 3 T MRI horizontal scanners closed bore or wide bore with a circularly polarised (CP) RF field for a maximum active scan time of 60 minutes.

WARNING

MR scans at 3 T must be performed in quadrature mode or CP mode for the radio frequency (RF) transmit coil. Using a multichannel mode may result in localised heating above safe levels.

A patient with one or two of these devices can be safely scanned in an MR system meeting conditions on the following pages.

Consider the following prior to scanning:

- Transmit/receive head coils and whole body coils may be safely used within the recommended SAR limits. Refer to the MRI safety information and recommended SAR limit tables in the following pages.
- Local cylindrical transmit/receive coils may be safely used, without SAR restriction, provided that the distance between the entire implant and the end of the local RF coil is at least equal to the radius of the local RF coil.
- It is safe to use local cylindrical RF receive only coils with implants during MRI scanning, provided SAR limits for the transmit coil have not been exceeded.
- Local planar (flat linearly polarized) receive only RF coils should be kept more than 10 cm away from the implant.

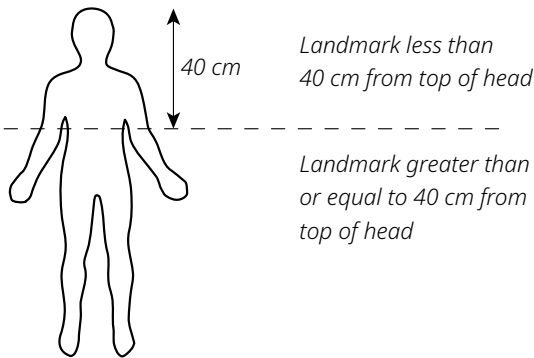


Fig.3: Landmark locations

OSI200 Implant and 1.5 T scans

- Remove the sound processor before entering the MRI scan room. The sound processor is MR Unsafe.
- Use the Cochlear Osia MRI Kit for MR scans at 1.5 T with the implant magnet in place. For instructions, refer to the Cochlear™ Osia® MRI Kit User Guide.
- Static magnetic field of 1.5 T.
- Maximum spatial field gradient of 2000 gauss/cm (20 T/m).
- When using a transmit/receive head coil, a maximum MR system reported, head averaged specific absorption rate (SAR) of 3.2 W/kg.
- When using a transmit body coil, a maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2 W/kg.

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant when imaged with a gradient echo pulse sequence scan in the axial plane is as follows:

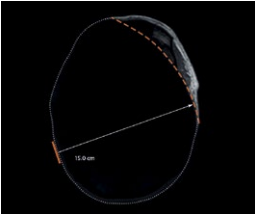
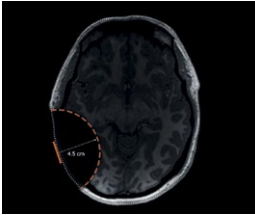
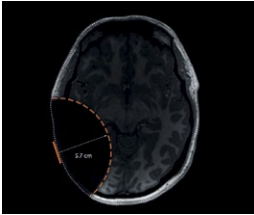
With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
		
15.0 cm (5.9 in)	4.5 cm (1.8 in)	5.7 cm (2.2 in)

Table 2: Maximum image artefact from centre at 1.5 T (gradient echo sequence). The image artefact may extend further in the coronal or sagittal plane.

NOTE

The image artefact results are based on worst-case scenarios showing maximum artefact extension. The optimisation of scan parameters can be used to minimise the extent of the artefact.

For bilateral OSI200 Implant recipients, the image artefacts as shown above are mirrored on the opposite side of the head for each implant. There may be some extension of the artefact between the implants.

OSI200 Implant and 3 T scans

- Surgically remove the implant magnet before MR scans at 3 T.
See **OSI200 Implant Physician's Guide** for additional information.
- Remove the sound processor before entering the MRI scan room. The sound processor is MR Unsafe.
- Static magnetic field of 3 T with the implant magnet surgically removed.
- Maximum spatial field gradient of 2000 gauss/cm (20 T/m).
- When using a transmit/receive head coil, a maximum MR system reported, head averaged specific absorption rate (SAR) of 3.2 W/kg.
- When using a transmit body coil, a maximum MR system reported, whole body averaged specific absorption rate (SAR) of 2 W/kg.
- Scans must be performed in circular polarization mode.

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant when imaged with a gradient echo pulse sequence scan in the axial plane is as follows:

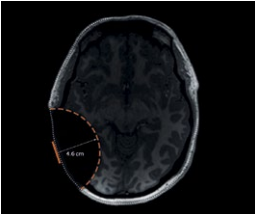
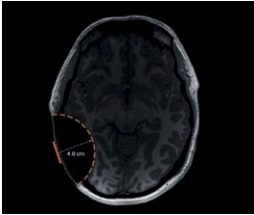
With non-magnetic plug	Implant magnet removed
	
4.6 cm (1.8 in)	4.0 cm (1.6 in)

Table 3: Maximum image artefact from centre at 3 T (gradient echo sequence). The image artefact may extend further in the coronal or sagittal plane.

Image interference and artefacts

The Cochlear Osia OSI200 Implant will create shadowing on the MR image near the implant, resulting in a loss of diagnostic information.

If inspecting near the implant, consider removing the implant magnet as MR image quality may be compromised with it in place.

If the implant magnet needs to be removed, refer the patient to an appropriate physician to arrange for the magnet to be removed before the MR scan.

The further optimisation of scan parameters can be used to minimise the extent of the artefact.

The image artefact extends from the centre of the implant. The Metal Artefact Reduction Sequence (MARS) parameters detailed in the tables below were used to produce the artefact sizes detailed in the following pages.

Parameter	Gradient echo	SEMAC-VAT	MARS	Optimized MARS
Scanning sequence	Gradient echo	Spin echo	Spin echo	Spin echo
Slice selection	Axial	Axial	Axial	Axial
Slice thickness	4 mm	3.5 mm	3 mm	3 mm
Repetition time	100 ms	3810 ms	4056 ms	3190 ms
Echo time	15 ms	80 ms	80 ms	80 ms
Echo train length	1	22	15	15
Pixel bandwidth	35 Hz/pixel	436 Hz/pixel	435 Hz/pixel	859 Hz/pixel
Acquisition matrix	256x256	343x310	499x451	499x442
Flip angle	30°	90°	90°	90°
SAR	0.02 W/kg	1.67 W/kg	1.20 W/kg	1.90 W/kg
dB/dt	6.86 T/s	71.64 T/s	88.40 T/s	93.38 T/s
B1rms	0.39 μ T	3.43 μ T	2.91 μ T	3.66 μ T
Duration	544 s (9min04s)	720 s (12 min)	709 s (11min49s)	344 s (5min44s)

Table 4: Scan parameters for scanning in a 1.5 T scanner

The following image artefact results are based on maximum artefact extension from the centre of the implant when scanned at 1.5 T using different Metal Artefact Reduction Sequence (MARS) sequences.

For bilateral implant recipients, the image artefacts as shown below are mirrored on the opposite side of the head for each implant. There may be some extension of the artefact between the implants.

OSI200 Implant and 1.5 T scans with SEMAC-VAT sequence

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant when imaged with a SEMAC-VAT sequence scan in the axial plane is:

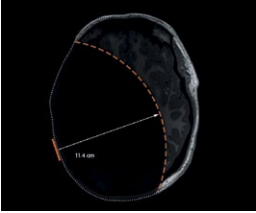
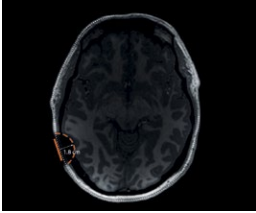
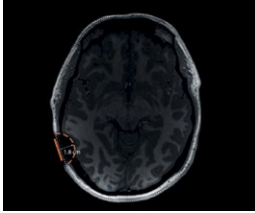
With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
		
11.4 cm (4.5 in)	1.8 cm (0.7 in)	1.8 cm (0.7 in)

Table 5: Maximum image artefact from centre at 1.5 T (SEMAC-VAT sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with a SEMAC-VAT sequence scan in the coronal plane is:

With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
7.5 cm (3.0 in)	2.0 cm (0.8 in)	2.0 cm (0.8 in)

Table 6: Maximum image artefact from centre at 1.5 T (SEMAC-VAT sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with a SEMAC-VAT sequence scan in the sagittal plane is:

With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
9.0 cm (3.5 in)	2.8 cm (1.1 in)	2.5 cm (1.0 in)

Table 7: Maximum image artefact from centre at 1.5 T (SEMAC-VAT sequence).

OSI200 Implant and 1.5 T scans with MARS sequence

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant when imaged with a MARS sequence scan in the axial plane is:

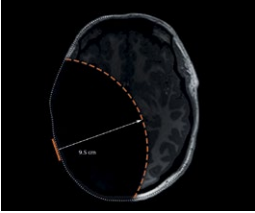
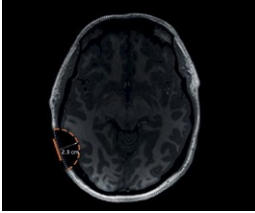
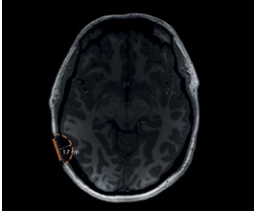
With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
		
9.5 cm (3.7 in)	2.3 cm (0.9 in)	1.7 cm (0.7 in)

Table 8: Maximum image artefact from centre at 1.5 T (MARS sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with a MARS sequence scan in the coronal plane is:

With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
6.6 cm (2.6 in)	2.3 cm (0.9 in)	1.9 cm (0.7 in)

Table 9: Maximum image artefact from centre at 1.5 T (MARS sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with a MARS sequence scan in the sagittal plane is:

With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
7.1 cm (2.8 in)	3.0 cm (1.2 in)	2.8 cm (1.1 in)

Table 10: Maximum image artefact from centre at 1.5 T (MARS sequence).

OSI200 Implant and 1.5 T scans with Optimized MARS sequence

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant when imaged with an Optimized MARS sequence scan in the axial plane is:

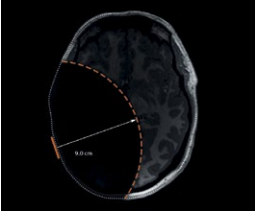
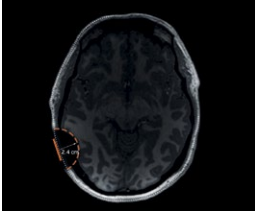
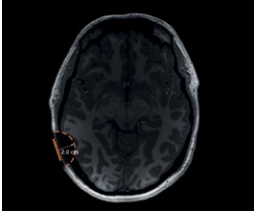
With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
		
9.0 cm (3.5 in)	2.4 cm (0.9 in)	2.0 cm (0.8 in)

Table 11: Maximum image artefact from centre at 1.5 T (Optimized MARS sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with an Optimized MARS sequence scan in the coronal plane is:

With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
8.2 cm (3.2 in)	1.9 cm (0.8 in)	1.7 cm (0.7 in)

Table 12: Maximum image artefact from centre at 1.5 T (Optimized MARS sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with an Optimized MARS sequence scan in the sagittal plane is:

With implant magnet + magnetic splint	With non-magnetic plug	Implant magnet removed
6.8 cm (2.7 in)	2.4 cm (0.9 in)	3.2 cm (1.3 in)

Table 13: Maximum image artefact from centre at 1.5 T (Optimized MARS sequence).

Parameter	Gradient echo	SEMAC-VAT	MARS
Scanning sequence	Gradient echo	Spin echo	Spin echo
Slice selection	Axial	Axial	Axial
Slice thickness	4 mm	3.5 mm	3 mm
Repetition time	100 ms	3197 ms	4809 ms
Echo time	9 ms	80 ms	80 ms
Echo train length	1	14	12
Pixel bandwidth	217 Hz/pixel	1244 Hz/pixel	1029 Hz/pixel
Acquisition matrix	256x256	307x277	300x268
Flip angle	80°	90°	90°
SAR	0.89 W/kg	1.88 W/kg	0.98 W/kg
dB/dt	16.20 T/s	58.31 T/s	53.21 T/s
B1rms	1.33 µT	1.93 µT	1.40 µT
Duration	182 s (3min02s)	409 s (6min49s)	289 s (4min49s)

Table 14: Scan parameters for scanning in a 3 T scanner

The following image artefact results are based on maximum artefact extension from the centre of the implant when scanned at 3 T using different Metal Artefact Reduction Sequence (MARS) sequences.

For bilateral implant recipients, the image artefacts as shown below are mirrored on the opposite side of the head for each implant. There may be some extension of the artefact between the implants.

OSI200 Implant and 3 T scans with SEMAC-VAT sequence

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant when imaged with a SEMAC-VAT sequence scan in the axial plane is:

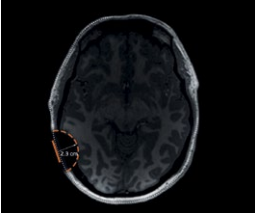
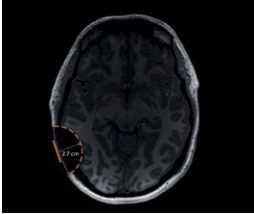
With non-magnetic plug	Implant magnet removed
	
2.3 cm (0.9 in)	2.7 cm (1.1 in)

Table 15: Maximum image artefact from centre at 3 T (SEMAC-VAT sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with an SEMAC-VAT sequence scan in the coronal plane is:

With non-magnetic plug	Implant magnet removed
2.4 cm (0.9 in)	2.4 cm (0.9 in)

Table 16: Maximum image artefact from centre at 3 T (SEMAC-VAT sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with an SEMAC-VAT sequence scan in the sagittal plane is:

With non-magnetic plug	Implant magnet removed
3.3 cm (1.3 in)	3.1 cm (1.2 in)

Table 17: Maximum image artefact from centre at 3 T (SEMAC-VAT sequence).

OSI200 Implant and 3 T scans with MARS sequence

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant when imaged with a MARS sequence scan in the axial plane is:

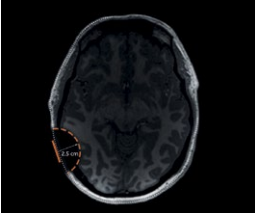
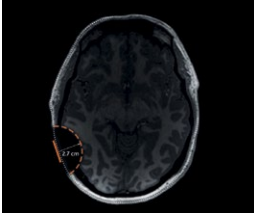
With non-magnetic plug	Implant magnet removed
	
2.5 cm (1.0 in)	2.7 cm (1.1 in)

Table 18: Maximum image artefact from centre at 3 T (MARS sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with a MARS sequence scan in the coronal plane is:

With non-magnetic plug	Implant magnet removed
2.9 cm (1.1 in)	2.6 cm (1.0 in)

Table 19: Maximum image artefact from centre at 3 T (MARS sequence).

In non-clinical testing, the maximum image artefact caused by the OSI200 Implant with a MARS sequence scan in the sagittal plane is:

With non-magnetic plug	Implant magnet removed
3.3 cm (1.3 in)	3.7 cm (1.5 in)

Table 20: Maximum image artefact from centre at 3 T (MARS sequence).

Considerations after an MRI examination

With the implant magnet in place

After the patient leaves the MRI room, remove the MRI kit bandage and splint. Ask the patient to place the sound processor on their head and turn it on. Confirm that the placement of the sound processor is correct and that there is no discomfort and sound is perceived as normal.

If there is discomfort or a change in sound perception, or problems with the placement of the sound processor, ask the patient to seek assistance from their implant clinician as soon as possible.

With the implant magnet removed

See *"Considerations for implant magnet removal"* on [page 7](#).

Disposal

The Cochlear Osia MRI Kit can be disposed of as normal hospital or household waste, or in accordance with local regulations.

The MRI Kit is for single-use only.

Notes

Notes



P2123908

AU Cochlear Ltd (ABN 96 002 618 073)
1 University Avenue, Macquarie University, NSW 2109,
Australia
Tel: +61 2 9428 6555 Fax: +61 2 9428 6352

ECREP DE Cochlear Deutschland GmbH & Co. KG
Mailänder Straße 4 a, 30539 Hannover, Germany
Tel: +49 511 542 770 Fax: +49 511 542 7720

CHIREP CH Cochlear AG
Peter Merian-Weg 4, 4052 Basel, Switzerland
Tel: +41 61 205 8204 Fax: +41 61 205 8205

US Cochlear Americas
10350 Park Meadows Drive, Lone Tree, CO 80124, USA
Tel: +1 303 790 9010

CA Cochlear Canada Inc
2500-120 Adelaide Street West, Toronto, ON M5H 1T1,
Canada
Tel: +1 (800) 483 3123 Fax: +1 416 972 5083

GB UK Responsible Person: Cochlear Europe Ltd
6 Dashwood Lang Road, Bourne Business Park, Addlestone,
Surrey KT15 2HJ, United Kingdom
Tel: +44 1932 26 3400 Fax: +44 1932 26 3426

BE Cochlear Benelux NV
Schaliënhoevdreef 20 i, B-2800 Mechelen, Belgium
Tel: +32 15 79 55 11 Fax: +32 15 79 55 70

FR Cochlear France S.A.S.
135 Route de Saint-Simon, 31035 Toulouse, France
Tel: +33 5 34 63 85 85 (International) or 0805 200 016
(National)
Fax: +33 5 34 63 85 80

IT Cochlear Italia S.r.l.
Via Trattati Comunitari Europei 1957-2007 n.17,
40127 Bologna (BO), Italy
Tel: +39 051 601 53 11 Fax: +39 051 39 20 62

SE Cochlear Nordic AB
Konstruktionsvägen 14, 435 33 Mölnlycke, Sweden
Tel: +46 31 335 14 61 Fax: +46 31 335 14 60

TR Cochlear Tibbi Cihazlar ve Sağlık Hizmetleri Ltd. Şti.
Küçükbakkalköy Mah, Defne Sok, Büyükhanlı Plaza No:3 Kat:3
Daire: 9-10-11-12, 34750, Ataşehir, İstanbul, Türkiye
Tel: +90 216 538 5900 Fax: +90 216 538 5919

HK Cochlear (HK) Limited
Room 1404-1406, 14/F, Leighton Centre, 77 Leighton Road,
Causeway Bay, Hong Kong
Tel: +852 2530 5773 Fax: +852 2530 5183

KR Cochlear Korea Ltd
2nd Floor, Yongsan Centreville Asterium, 25,
Hangang-daero 30 gil, Yongsan-gu, Seoul, Korea (04386)
Tel: +82 2 533 4450 Fax: +82 2 533 8408

CN Cochlear Medical Device (Beijing) Co., Ltd
Unit 2608-2617, 26th Floor, No.9 Building, No.91 Jianguo
Road,
Chaoyang District, Beijing 100022, P.R. China
Tel: +86 10 5909 7800 Fax: +86 10 5909 7900

IN Cochlear Medical Device Company India Pvt. Ltd.
Ground Floor, Platina Building, Plot No C-59, G-Block,
Bandra Kurla Complex, Bandra (E), Mumbai - 400 051, India
Tel: +91 22 6112 1111 Fax: +91 22 6112 1100

JP 株式会社日本コクレア(Nihon Cochlear Co Ltd)
〒113-0033 東京都文京区本郷2-3-7 お茶の水元町ビル
Tel: +81 3 3817 0241 Fax: +81 3 3817 0245

AE Cochlear Middle East FZ-LLC
Dubai Healthcare City, Al Razi Building 64, Block A, Ground
Floor, Offices IR1 and IR2, Dubai, United Arab Emirates
Tel: +971 4 818 4400 Fax: +971 4 361 8925

PA Cochlear Latinoamérica S.A.
International Business Park, Building 3835, Office 403,
Panama Pacifico, Panama
Tel: +507 830 6220 Fax: +507 830 6218

NZ Cochlear NZ Limited
Level 4, Takapuna Towers, 19-21 Como St, Takapuna,
Auckland 0622, New Zealand
Tel: + 64 9 914 1983 Fax: 0800 886 036

www.cochlear.com

This material is intended for health professionals. If you are a consumer, please seek advice from your health professional about treatments for hearing loss. Outcomes may vary, and your health professional will advise you about the factors which could affect your outcome. Always read the instructions for use. Not all products are available in all countries. Please contact your local Cochlear representative for product information.

Cochlear, 科利耳, コクレア, 코클리어, Hear now. And always, Osia, SmartSound, the elliptical logo, and marks bearing an ® or ™ symbol, are either trademarks or registered trademarks of the Cochlear group of companies (unless otherwise noted).

© Cochlear Limited 2023. All rights reserved. 2023-09

P1752884 D1752875-V3

CE
0123



D1752875-V3